	Test Case Name	Approve Sample Acknowledgement
	Test Case Number	LIMS-SM-014
	Actual Result	Pass

1. Enter the username, password, Select Module and click on login button

SSNo : PSS-LIMS-SS-001





2. Click on Sample Acknowledgement

SSNo : PSS-LIMS-SS-002

PHARMA SOFT

STABILITY MANAGEMENT



LIMS Approver

Sample Acknowledgment Approval

Activity :
 . Select The Record and Click On view.
 . Enter Comments and Click on Approve / Revert.
 . Enter E-Signature Password and Click on Submit.

S.No	Protocol Number	Product Name	Specification Number	AR Number	Storage Condition	Stability Condition	Packing type	Batch Number	Start Date	Done By/On	Action
1	PC0006	PARA65	SPSSPARA65	ASS220042	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	Blister	Product Batch No	16/08/2022	LIMS User /16/08/2022 19:15:55	View
2	PC0007	PARA65	SPSSPARA65	ASS220050	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	Blister	Product Batch No 001	17/08/2022	LIMS Approver /17/08/2022 14:53:43	View
3	PC0018	PARA65	SPPARACH001	ASS220080	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	Blister	Batch001	01/09/2022	LIMS User /01/09/2022 15:17:08	View

<< < 1 > >>


Go to page:

Row count:

Showing 1-3 of 3

LIMS 2.0

All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Approve Sample Acknowledgement
	Test Case Number	LIMS-SM-014
	Actual Result	Pass

3. Select a Record And Click on View

SSNo : PSS-LIMS-SS-003

PHARMA SOFT
STABILITY MANAGEMENT

LIMS Approver

Sample Acknowledgement

Activity :
. Select The
. Enter Com
. Enter E-Sig

List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

Protocol Number: PC0018
Product Name: PARA65

Specification Number: SPPARACH001
AR Number: ASS220080

Batch Number: Batch001
Packing Type: Blister

Storage Conditions: 40°C ± 2°C / 75%RH ± 5%RH
Interval: MONTH

Qty Withdraw: 2 UOM: Gm
Total Qty Remained: 43 UOM: Gm

Start Date: 01/09/2022

Used For: Sampling

Comments*: Maximum allowed characters should be 1500

Approve Revert

Action

View

View

View

Showing 1-3 of 3

4. Enter Comments

SSNo : PSS-LIMS-SS-004

PHARMA SOFT
STABILITY MANAGEMENT

LIMS Approver

Sample Acknowledgement

Activity :
. Select The
. Enter Com
. Enter E-Sig

List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

Protocol Number: PC0018
Product Name: PARA65

Specification Number: SPPARACH001
AR Number: ASS220080

Batch Number: Batch001
Packing Type: Blister

Storage Conditions: 40°C ± 2°C / 75%RH ± 5%RH
Interval: MONTH

Qty Withdraw: 2 UOM: Gm
Total Qty Remained: 43 UOM: Gm

Start Date: 01/09/2022

Used For: Sampling

Comments*: Approved

Approve Revert


Action

View

View


View

Showing 1-3 of 3

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Approve Sample Acknowledgement
	Test Case Number	LIMS-SM-014
	Actual Result	Pass

5. Click on Approve

SSNo : PSS-LIMS-SS-005

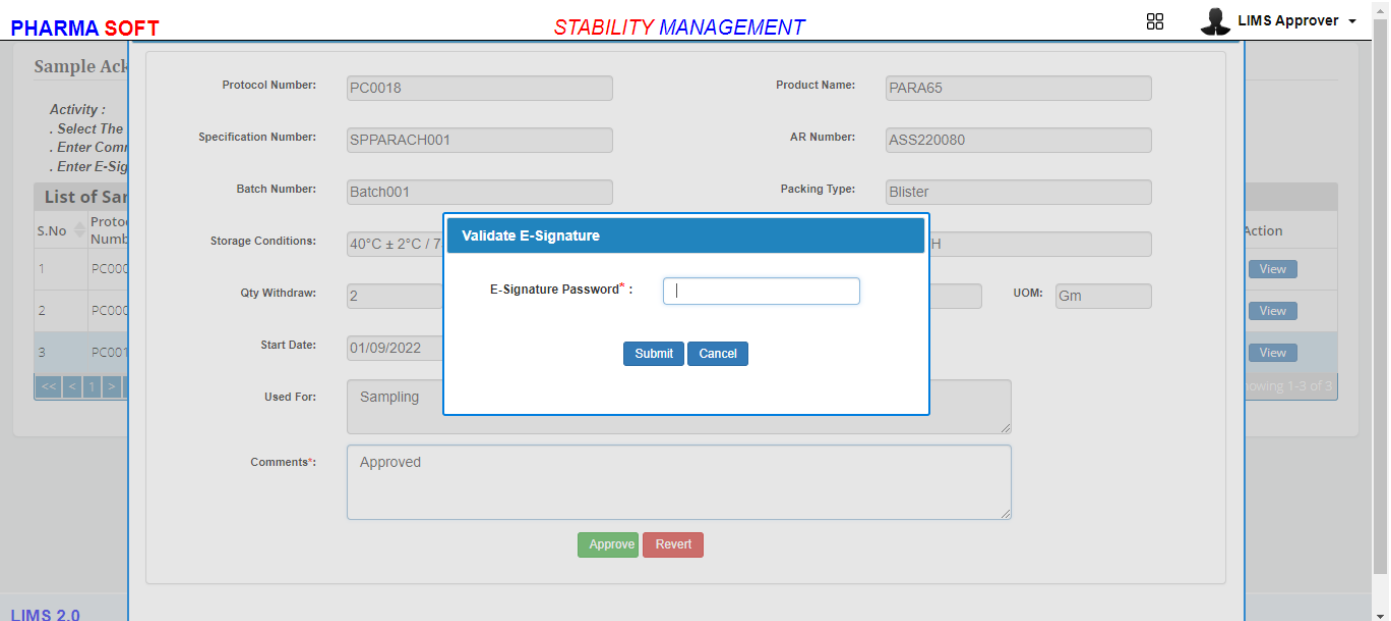

PHARMA SOFT **STABILITY MANAGEMENT** LIMS Approver

Sample Acknowledgement

Activity :
 . Select The
 . Enter Com
 . Enter E-Sig

List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

Protocol Number: PC0018 Product Name: PARA65

Specification Number: SPPARACH001 AR Number: ASS220080

Batch Number: Batch001 Packing Type: Blister

Storage Conditions: 40°C ± 2°C / 7 H

Qty Withdraw: 2 UOM: Gm

Start Date: 01/09/2022

Used For: Sampling

Comments*: Approved

Validate E-Signature

E-Signature Password* :


Submit Cancel

Approve Revert

LIMS 2.0

6. Enter E-Signature Password

SSNo : PSS-LIMS-SS-006

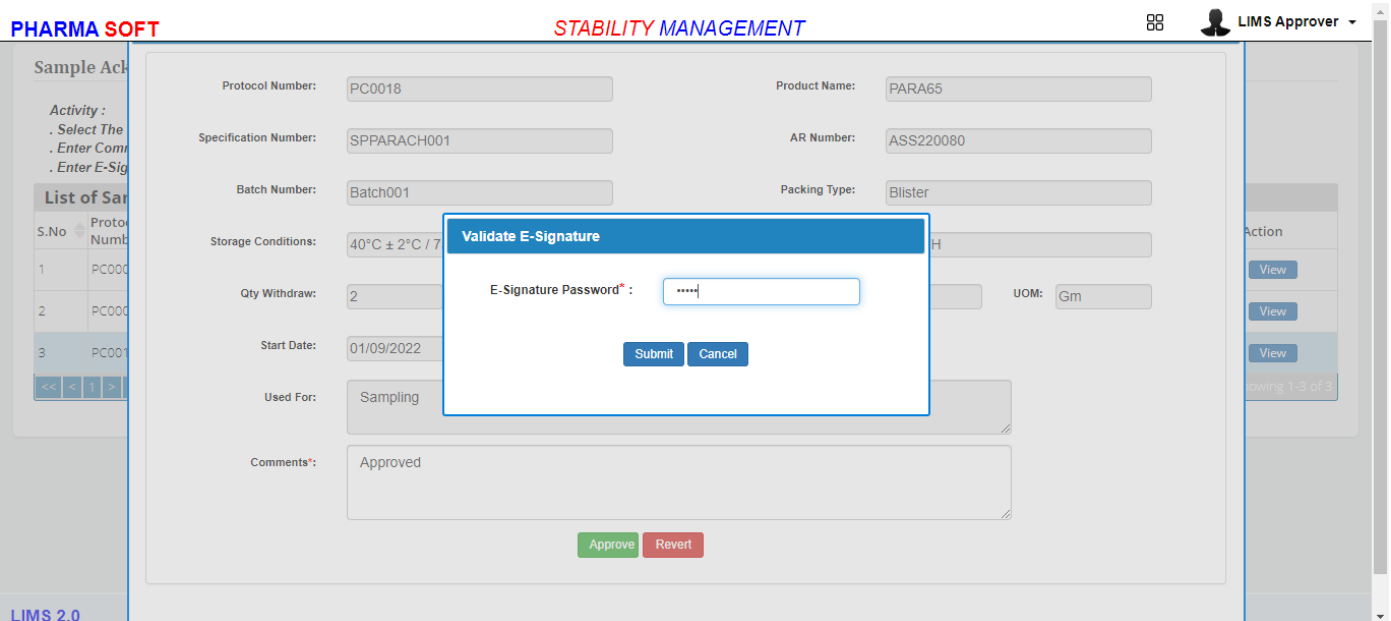

PHARMA SOFT **STABILITY MANAGEMENT** LIMS Approver

Sample Acknowledgement

Activity :
 . Select The
 . Enter Com
 . Enter E-Sig

List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

Protocol Number: PC0018 Product Name: PARA65

Specification Number: SPPARACH001 AR Number: ASS220080

Batch Number: Batch001 Packing Type: Blister

Storage Conditions: 40°C ± 2°C / 7 H

Qty Withdraw: 2 UOM: Gm

Start Date: 01/09/2022

Used For: Sampling

Comments*: Approved


Validate E-Signature

E-Signature Password* :

Submit Cancel

Approve Revert

LIMS 2.0

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Approve Sample Acknowledgement
	Test Case Number	LIMS-SM-014
	Actual Result	Pass

7. Click on Submit

SSNo : PSS-LIMS-SS-007

PHARMA SOFT
STABILITY MANAGEMENT
☰
LIMS Approver ▾

Sample Acknowledgement

Activity :
 . Select The
 . Enter Com
 . Enter E-Sig

List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

Protocol Number: PC0018
Product Name: PARA65

Specification Number: SPPARACH001
AR Number: ASS220080

Batch Number: Batch001
Packing Type: Blister

Storage Conditions: 40°C ± 2°C / 75%RH ± 5%RH
Interval: MONTH

Qty Withdraw: 2 UOM: Gm
Total Qty Remained: 43 UOM: Gm

Start Date: 01/09/2022

Used For: Sampling

Comments*: Approved

Approve Revert

View
View
View

Showing 1-3 of 3

8. Click on OK button

SSNo : PSS-LIMS-SS-008

PHARMA SOFT
STABILITY MANAGEMENT
☰
LIMS Approver ▾

Sample Acknowledgement

Activity :
 . Select The
 . Enter Com
 . Enter E-Sig

List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

Protocol Number: PC0018
Product Name: PARA65

Specification Number: SPPARACH001
AR Number: ASS220080

Batch Number: Batch001
Packing Type: Blister

Storage Conditions: 40°C ± 2°C / 75%RH ± 5%RH
Interval: MONTH

Qty Withdraw: 2 UOM: Gm
Total Qty Remained: 43 UOM: Gm

Start Date: 01/09/2022


Used For: Sampling

Comments*: Approved

Approve Revert

View
View
View

Showing 1-3 of 3

	Test Case Name	Approve Sample Acknowledgement
	Test Case Number	LIMS-SM-014
	Actual Result	Pass


9.Click on User Name

SSNo : PSS-LIMS-SS-009

PHARMA SOFT SOL
STABILITY MANAGEMENT
LIMS Approver

Sample Acknowledgment Approval

Activity :
 . Select The Record and Click On view.
 . Enter Comments and Click on Approve / Revert.
 . Enter E-Signature Password and Click on Submit.



List of Sample Acknowledgments for Approval


S.No	Protocol Number	Product Name	Specification Number	AR Number	Storage Condition	Stability Condition	Packing type	Batch Number	Start Date	Done By/On	Log	Action
1	PC0006	PARA65	SPSSPARA65	ASS220042	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	Blister	Product Batch No	16/08/2022	LIMS User /16/08/2022 19:15:55		View
2	PC0007	PARA65	SPSSPARA65	ASS220050	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	Blister	Product Batch No 001	17/08/2022	LIMS Approver /17/08/2022 14:53:43		View

<< < 1 > >> Go to page: Row count: Showing 1-2 of 2

LIMS 2.0
All © rights Reserved PharmaSoftSol

10.Click on Logout

SSNo : PSS-LIMS-SS-010



LIMS

User Name

Password

Module