
	Test Case Name	Create Schedule Test
	Test Case Number	LIMS-SM-007
	Actual Result	Pass

1.Enter the username, password, Select Module and click on login button

SSNo : PSS-LIMS-SS-001

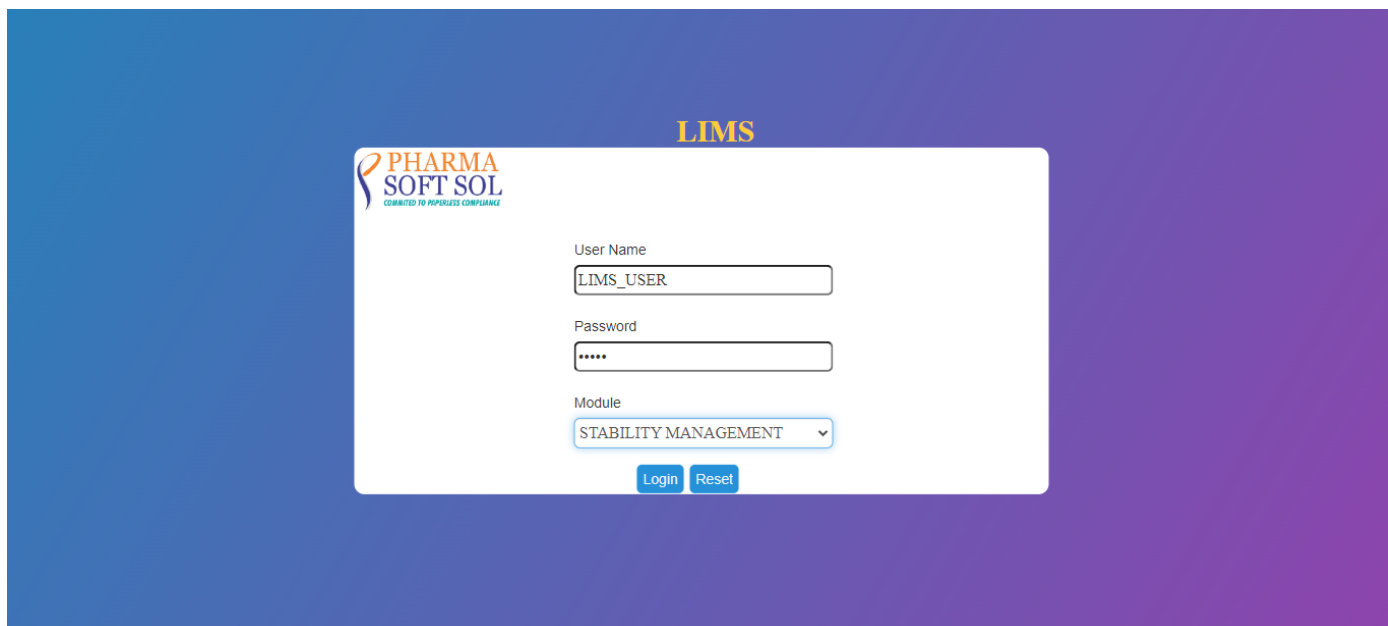


The screenshot shows the LIMS login interface. At the top center, the word "LIMS" is displayed in yellow. Below it is the PharmaSoftSol logo. The login form contains three input fields: "User Name" with the value "LIMS_USER", "Password" with masked characters "*****", and "Module" with a dropdown menu set to "STABILITY MANAGEMENT". At the bottom of the form are two buttons: "Login" and "Reset".

2.Click on Schedule Test

SSNo : PSS-LIMS-SS-002



The screenshot displays the "Schedule Test Selection" screen within the "STABILITY MANAGEMENT" module. The page header includes the PharmaSoftSol logo, the module name "STABILITY MANAGEMENT", and a user profile icon labeled "LIMS USER". The main content area has three tabs: "Action", "Schedule Intervals", and "Schedule Test Selection", with the third tab being active. Below the tabs, an "Activity" section provides instructions: ". Select Create, Select The Record and Click on Next." and ". If Rejected is Required, Then Select Rejected." There are two radio buttons: "Schedule Test Selection" (unselected) and "Reverted Schedule Test Selection" (selected). Below the radio buttons is a table titled "List of Reverted Schedules". The table has columns: "S.No.", "Protocol Number", "Product Name", "Generic Name", "Specification Name", "Reverted By", "Reverted On", and "Comments". The table body is empty, with a message "No data available!". At the bottom left of the table, there is a "Row count:" label and a dropdown menu set to "10". The footer of the page shows "LIMS 2.0" on the left and "All © rights Reserved PharmaSoftSol" on the right.

	Test Case Name	Create Schedule Test
	Test Case Number	LIMS-SM-007
	Actual Result	Pass

3.Click on Create

SSNo : PSS-LIMS-SS-003

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Schedule Test Selection

Action
Schedule Intervals
Schedule Test Selection

Activity :
 . Select Create, Select The Record and Click on Next.
 . If Rejected is Required, Then Select Rejected.

☒ Schedule Test Selection
 ☐ Reverted Schedule Test Selection

Schedule Create

S.No.	Protocol Number	Product	Specification Name	Generic Name
No data available!				

Go to page: 1 Row count: 10

IMS 2.0
 All © rights Reserved PharmaSoftSol

4.Select a Record

SSNo : PSS-LIMS-SS-004

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Schedule Test Selection

Action
Schedule Intervals
Schedule Test Selection

Activity :
 . Select Create, Select The Record and Click on Next.
 . If Rejected is Required, Then Select Rejected.


☒ Schedule Test Selection
 ☐ Reverted Schedule Test Selection

Schedule Create

S.No.	Protocol Number	Product	Specification Name	Generic Name
1	PS0011	Aspirin	SS Aspirin Spec	Aspirin and codeine
2	PS0023	Aspirin	SS Aspirin Spec	Aspirin and codeine

<< < 1 > >> Go to page: 1 Row count: 10

IMS 2.0
 All © rights Reserved PharmaSoftSol

	Test Case Name	Create Schedule Test
	Test Case Number	LIMS-SM-007
	Actual Result	Pass

5. Click on Next

SSNo : PSS-LIMS-SS-005

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Schedule Test Selection

Action
Schedule Intervals
Schedule Test Selection

☐ Reverted Schedule Test Selection

Specification Name	Generic Name
SS Aspirin Spec	Aspirin and codeine
SS Aspirin Spec	Aspirin and codeine

Showing 1-2 of 2

. Select The Record and Click On Next.

Schedule Test Selection				
<input type="checkbox"/> S.No	Packing Type	Product Batch No	Mfg.Date	Exp.Date
Go to page: 1 Row count: 10				

IMS 2.0
All © rights Reserved PharmaSoftSol

6. Select a Record

SSNo : PSS-LIMS-SS-006

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Schedule Test Selection


Action
Schedule Intervals
Schedule Test Selection

. Select The Record and Click On Next.

<input checked="" type="checkbox"/> S.No	Packing Type	Product Batch No	Mfg.Date	Exp.Date	Storage Conditions	Orientation	Stability Conditions	Schedule Interval Type
<input checked="" type="checkbox"/> 1	Box	44	02-11-2020	02-11-2021	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	MONTH

<< < 1 > >> Go to page: 1 Row count: 10 Showing 1-1 of 1

IMS 2.0
All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Schedule Test
	Test Case Number	LIMS-SM-007
	Actual Result	Pass

7. Click on Next

SSNo : PSS-LIMS-SS-007

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Schedule Test Selection

Action
Schedule Intervals
Schedule Test Selection

. Select The Test, Select Approval Life Cycle and Click On Finish.
. Enter E-Signature Password and Click On Submit.

Set1

Packing Type:

Storage Conditions:

Stability Condition:

Product Batch No:

Schedule Interval Type:

Orientation:

Interval and Test Details		
S.No	Test Name	1
1	Descriptive	<input checked="" type="checkbox"/>
2	Loss On Dry	<input checked="" type="checkbox"/>
S.No		1
1	Number of Syringes required at Each Point	<input type="text"/>

Approval Life Cycle*: View

IMS 2.0
All © rights Reserved PharmaSoftSol

8. Select Approval Cycle

SSNo : PSS-LIMS-SS-008

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Schedule Test Selection

Action
Schedule Intervals
Schedule Test Selection

. Select The Test, Select Approval Life Cycle and Click On Finish.
. Enter E-Signature Password and Click On Submit.

Set1

Packing Type:

Storage Conditions:

Stability Condition:

Product Batch No:


Schedule Interval Type:

Orientation:

Interval and Test Details		
S.No	Test Name	1
1	Descriptive	<input checked="" type="checkbox"/>
2	Loss On Dry	<input checked="" type="checkbox"/>
S.No		1
1	Number of Syringes required at Each Point	<input type="text" value="50"/>


Approval Life Cycle*: View

IMS 2.0
All © rights Reserved PharmaSoftSol

	Test Case Name	Create Schedule Test
	Test Case Number	LIMS-SM-007
	Actual Result	Pass

9. Click on Finish

SSNo : PSS-LIMS-SS-009

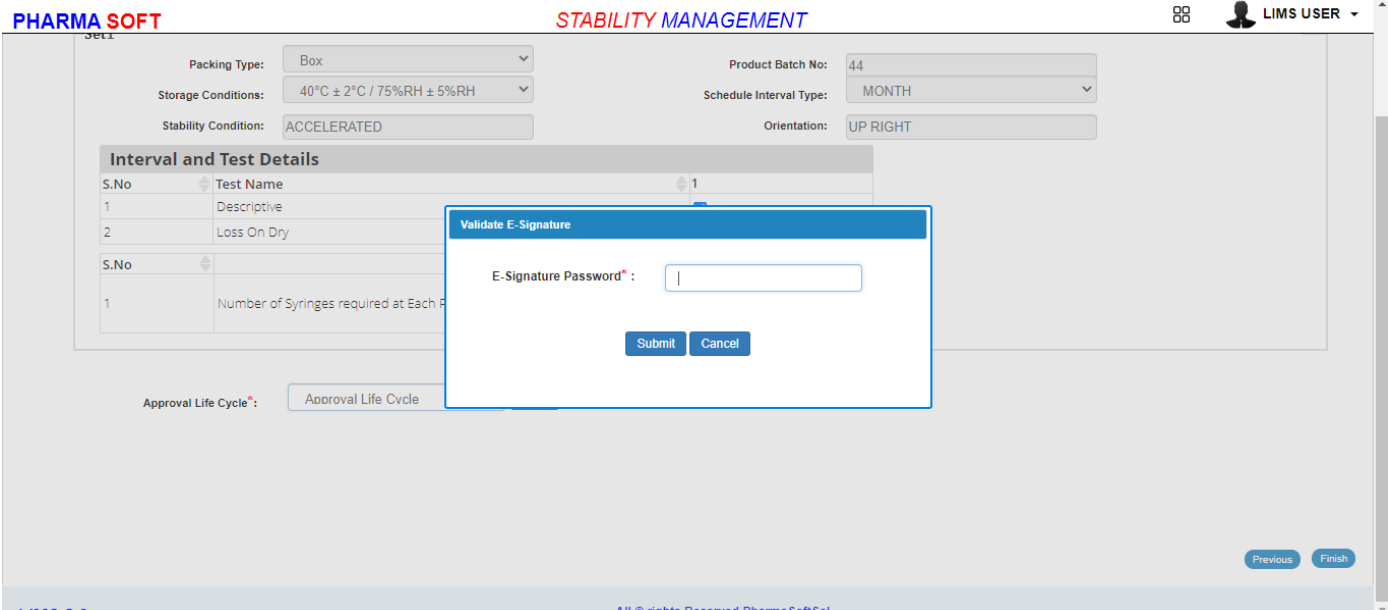

PHARMA SOFT **STABILITY MANAGEMENT** LIMS USER

Packing Type: Box Product Batch No: 44
Storage Conditions: 40°C ± 2°C / 75%RH ± 5%RH Schedule Interval Type: MONTH
Stability Condition: ACCELERATED Orientation: UP RIGHT

Interval and Test Details

S.No	Test Name
1	Descriptive
2	Loss On Dry

Validate E-Signature

E-Signature Password* :

Submit Cancel


Approval Life Cycle* : Approval Life Cycle

Previous Finish

LIMS 2.0 All © rights Reserved PharmaSoftSol

10. Enter E-Signature Password

SSNo : PSS-LIMS-SS-010

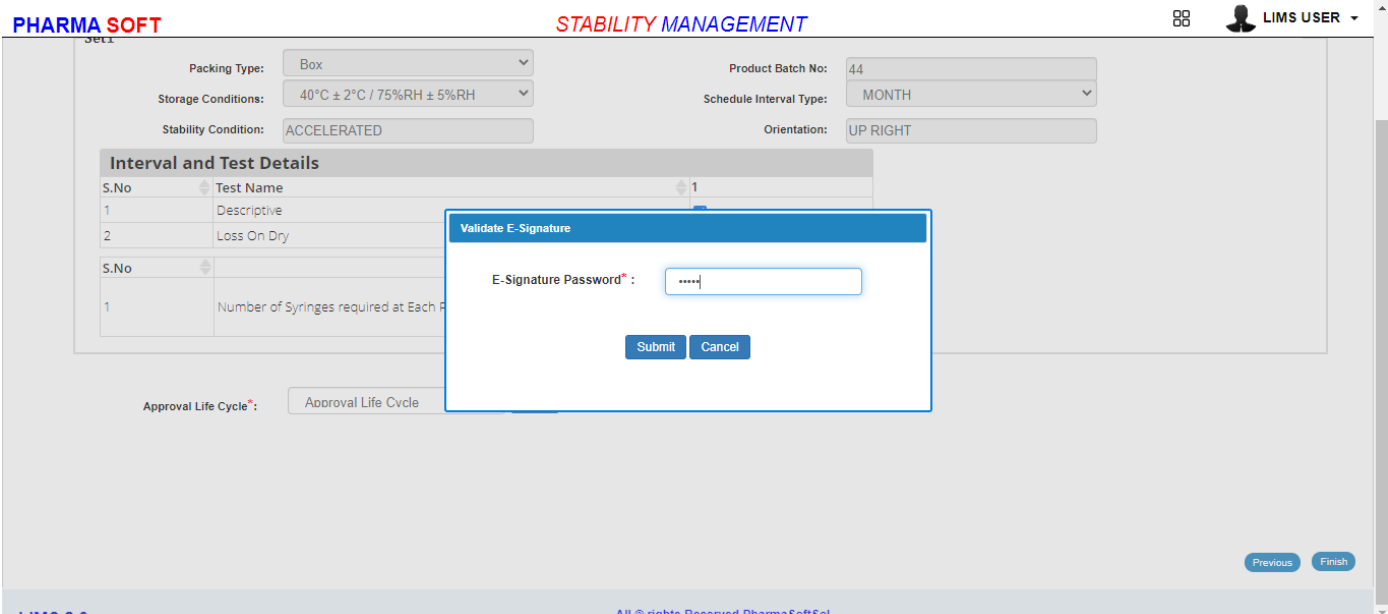

PHARMA SOFT **STABILITY MANAGEMENT** LIMS USER

Packing Type: Box Product Batch No: 44
Storage Conditions: 40°C ± 2°C / 75%RH ± 5%RH Schedule Interval Type: MONTH
Stability Condition: ACCELERATED Orientation: UP RIGHT

Interval and Test Details

S.No	Test Name
1	Descriptive
2	Loss On Dry

Validate E-Signature


E-Signature Password* :

Submit Cancel

Approval Life Cycle* : Approval Life Cycle

Previous Finish

LIMS 2.0 All © rights Reserved PharmaSoftSol

	Test Case Name	Create Schedule Test
	Test Case Number	LIMS-SM-007
	Actual Result	Pass

11.Click on OK button

SSNo : PSS-LIMS-SS-011

PHARMA SOFT **STABILITY MANAGEMENT** LIMS USER

Packing Type: Box Product Batch No: 44
Storage Conditions: 40°C ± 2°C / 75%RH ± 5%RH Schedule Interval Type: MONTH
Stability Condition: ACCELERATED Orientation: UP RIGHT

Interval and Test Details

S.No	Test Name
1	Descriptive
2	Loss On Dry

Success
Submitted successfully.
OK

Approval Life Cycle*: Approval Life Cycle view

Previous Finish

LIMS 2.0 All © rights Reserved PharmaSoftSol

12.Click on User Name

SSNo : PSS-LIMS-SS-012

PHARMA SOFT **STABILITY MANAGEMENT** LIMS USER

☐ Schedule Test Selection ☒ Reverted Schedule Test Selection

List of Reverted Schedules


S.No.	Protocol Number	Product Name	Generic Name	Specification Name	Reverted By	Reverted On	Comments
No data available!							

Row count: 10

Logout


Previous Next

LIMS 2.0 All © rights Reserved PharmaSoftSol


	Test Case Name	Create Schedule Test
	Test Case Number	LIMS-SM-007
	Actual Result	Pass

13. Click on Logout

SSNo : PSS-LIMS-SS-013

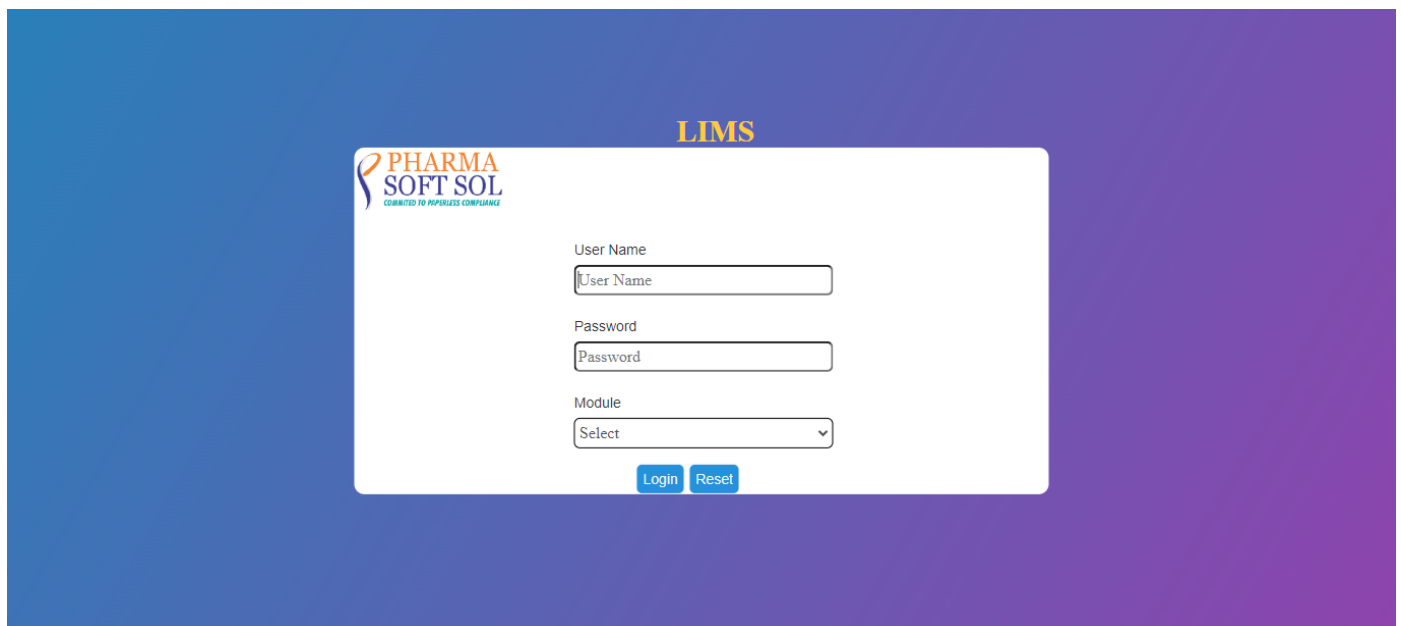

LIMS



User Name

Password

Module