
	Test Case Name	ReSubmit RejectedSampleAcknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass


1.Enter the username, password, Select Module and click on login button


SSNo : PSS-LIMS-SS-001

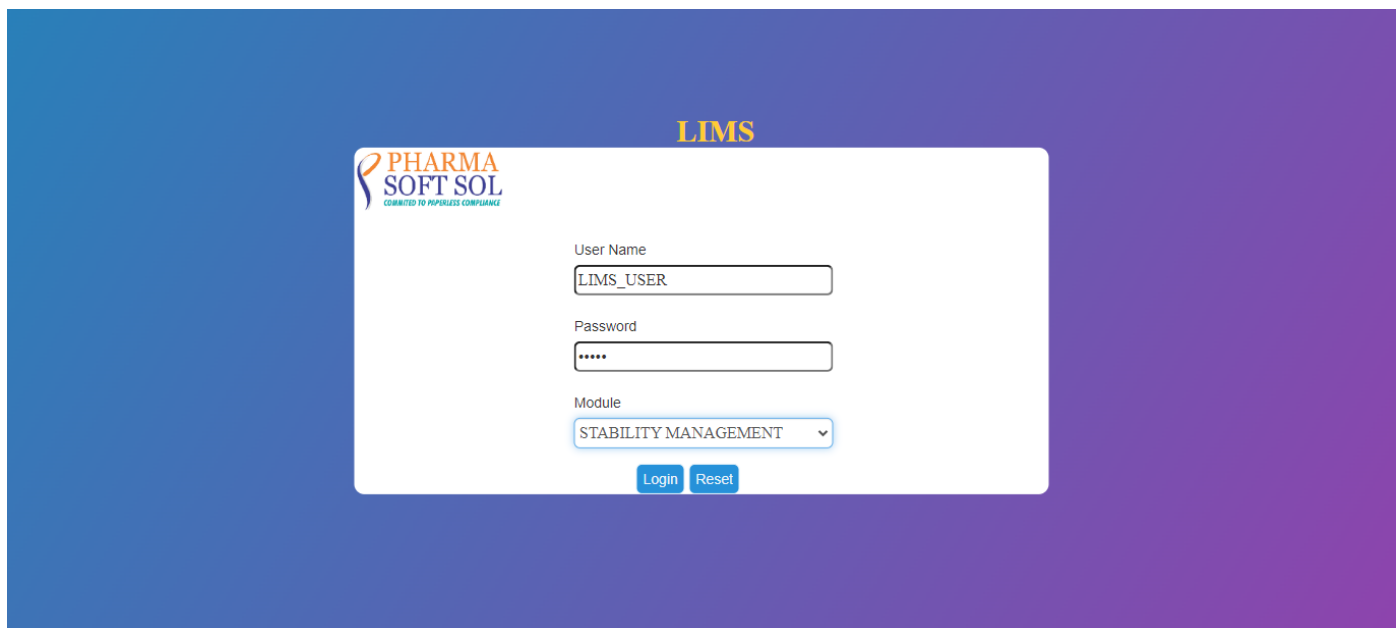

2.Click on Sample Acknowledgement

SSNo : PSS-LIMS-SS-002

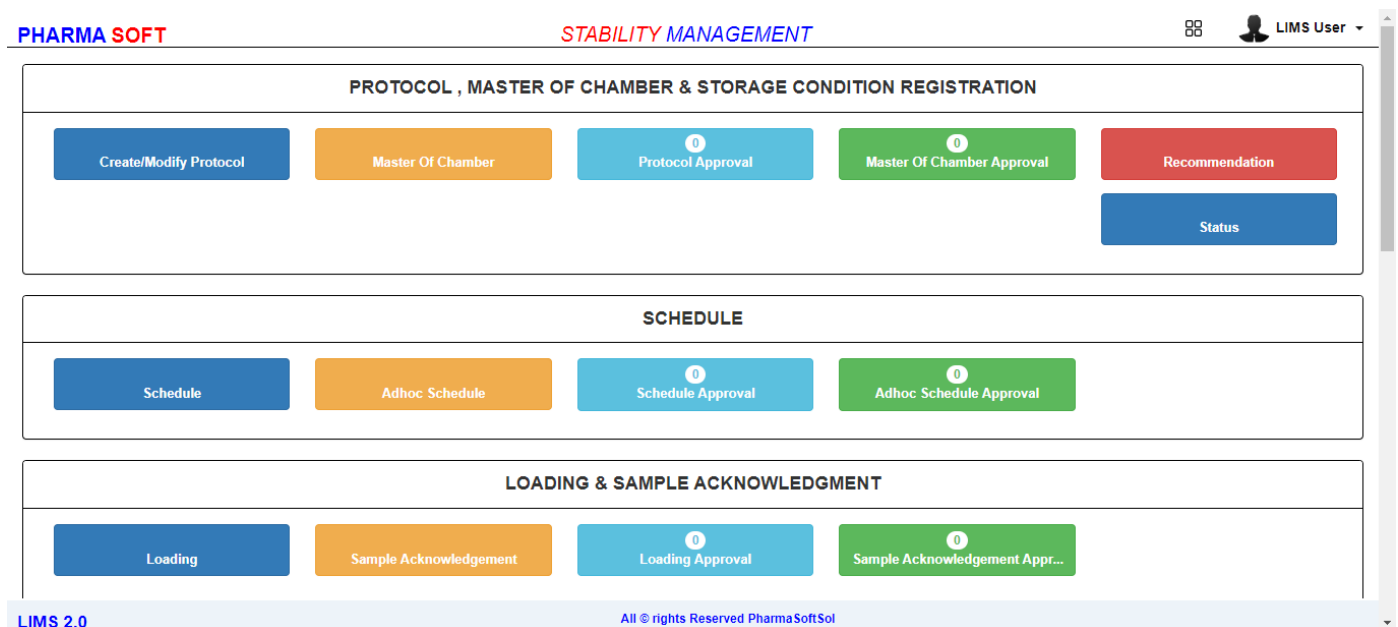

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	ReSubmit RejectedSampleAcknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

3.Select a Record

SSNo : PSS-LIMS-SS-003

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Sample Acknowledgment

Action
Sample Acknowledgment Details

. Select Sample Acknowledgment, Click On Search, Select The Record and Click on Next.
. If Reverted Sample Acknowledgments by Approver is Required, Then Select Reverted Sample Acknowledgments by Approver.

☐ Sample Acknowledgment
 ☒ Reverted Sample Acknowledgments by Approver

List of Reverted Sample Acknowledgements by Approver

S.No	Product Name	Specification Number	AR Number	Protocol Number	Storage Conditions	Orientation	Stability Conditions	Packing Type	Batch Number	Start Date	RevertedBy/On	Comments
1	PARA65	SPSSPARA65	ASS220051	PC0007	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	Blister	Product Batch No 001	17/08/2022	LIMS User 17/08/2022	view
2	PARA65	SPPARACH001	ASS220070	PC0015	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	Blister	44	31/08/2022	LIMS Approver 31/08/2022	view

<< < 1 > >>
Go to page: 1
Row count: 10
Showing 1-2 of 2

IMS 2.0
All © rights Reserved PharmaSoftSol

4.Click on Next

SSNo : PSS-LIMS-SS-004


PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Sample Acknowledgment

Action
Sample Acknowledgment Details

- Enter Quantity withdraw and select UOM from drop down.click on Select to select the Approver.
- If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish
- Enter E- Signature Password and Click On Finish.

IMS 2.0
All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	ReSubmit RejectedSampleAcknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

5. Click on Select Approver

SSNo : PSS-LIMS-SS-005

PHARMA SOFT

STABILITY MANAGEMENT

LIMS User

Sample Acknowledgement

Action

Sample Acknowledgement

• Enter Quantity

• If select Acknowledgement

• Enter E- Signature

Quantity

Withdrawal

Used For

Locations

organization

Plant-1

First Name:

Last Name:

Search

Clear

Users

First Name

Last Name

No data available!

Go to page: 1 Row count: 5

Select

Cancel

LIMS 2.0

All © rights Reserved PharmaSoftSol

6. Select Location

SSNo : PSS-LIMS-SS-006

PHARMA SOFT

STABILITY MANAGEMENT

LIMS User

Sample Acknowledgment

Action

Sample Acknowledgment

- Enter Quantity
- If select Acknowledgment
- Enter E- Signature

Quantity

Withdrawal

Used For

Locations

organization

Plant-1

First Name:

Last Name:

Search

Clear

Users

First Name

Last Name

No data available!

Go to page

1

Row count


5

Select

Cancel


LIMS 2.0

All © rights Reserved PharmaSoftSol

	Test Case Name	ReSubmit RejectedSampleAcknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

7.Select a Record

SSNo : PSS-LIMS-SS-007

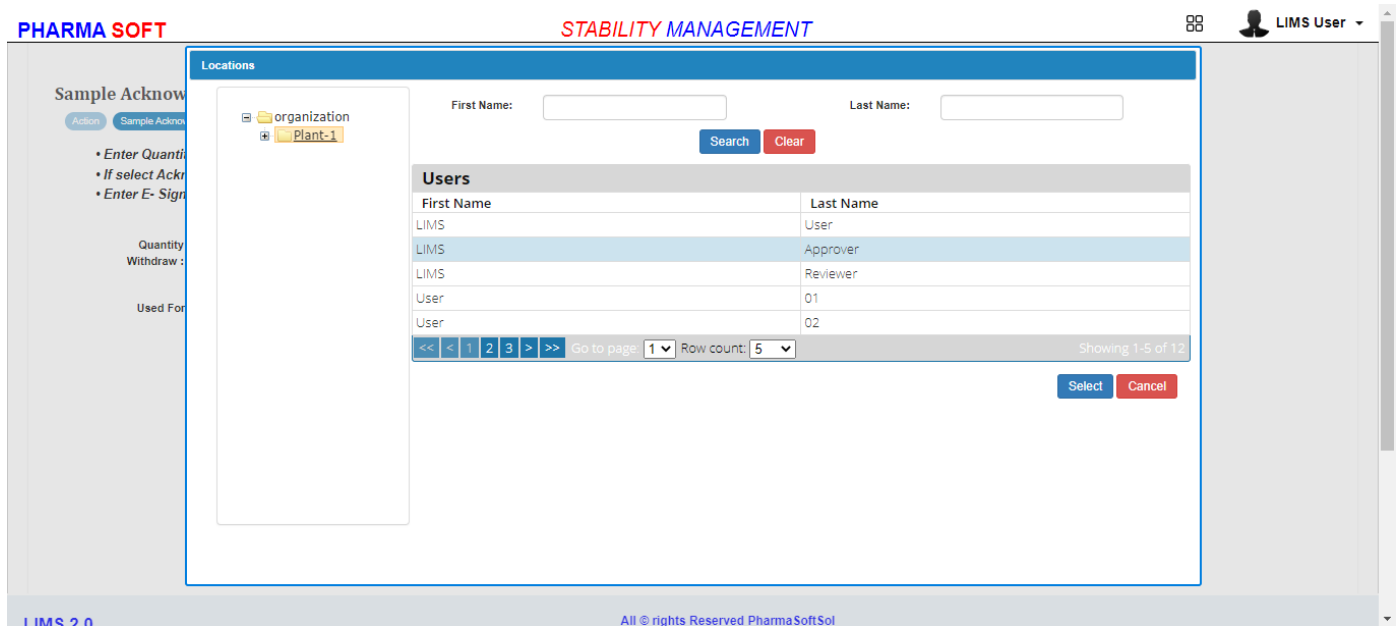


The screenshot shows the PHARMA SOFT STABILITY MANAGEMENT interface. The 'Locations' section on the left lists 'organization' and 'Plant-1'. The 'Users' section on the right displays a table of users with columns for First Name and Last Name. The table contains five rows of data, with the first row highlighted. Below the table, there are pagination controls showing 'Go to page: 1' and 'Row count: 5'. The 'Select' button is visible at the bottom right of the table.


First Name	Last Name
LIMS	User
LIMS	Approver
LIMS	Reviewer
User	01
User	02

8.Click on Select

SSNo : PSS-LIMS-SS-008

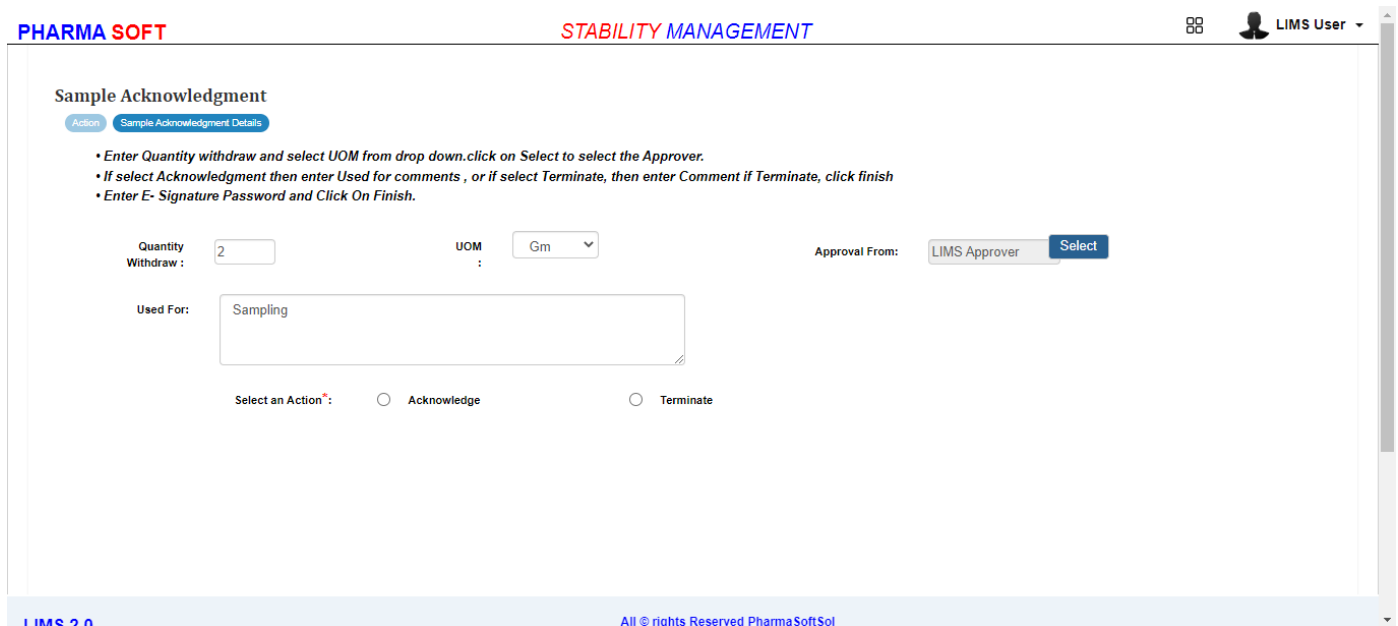

The screenshot shows the PHARMA SOFT STABILITY MANAGEMENT interface with the 'Sample Acknowledgment' form. The form includes fields for 'Quantity Withdraw' (set to 2), 'UOM' (set to Gm), and 'Approval From' (set to LIMS Approver). The 'Used For' field is set to 'Sampling'. Below the form, there are radio buttons for 'Acknowledge' and 'Terminate'. The 'Select' button is visible next to the 'Approval From' field.

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	ReSubmit RejectedSampleAcknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

9.Click on Acknowledgment Radio Button

SSNo : PSS-LIMS-SS-009

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Sample Acknowledgment

Action
Sample Acknowledgment Details

- Enter Quantity withdraw and select UOM from drop down.click on Select to select the Approver.
- If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish
- Enter E- Signature Password and Click On Finish.

Quantity Withdraw :
UOM :
Approval From:

Used For:

Select an Action*: ☒ Acknowledge ☐ Terminate

IMS 2.0
All © rights Reserved PharmaSoftSol

10.Click on Finish

SSNo : PSS-LIMS-SS-010

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Sample Acknowledgment

Action
Sample Acknowledgment Details

- Enter Quantity withdraw and select UOM from drop down.click on Select to select the Approver.
- If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish
- Enter E- Signature Password and Click On Finish.


Quantity Withdraw :
UOM :
Approval From:

Used For:

Select an Action*: ☒ Acknowledge ☐ Terminate


Validate E-Signature
E-Signature Password*:

IMS 2.0
All © rights Reserved PharmaSoftSol

	Test Case Name	ReSubmit RejectedSampleAcknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

11. Enter E-Signature Password


SSNo : PSS-LIMS-SS-011

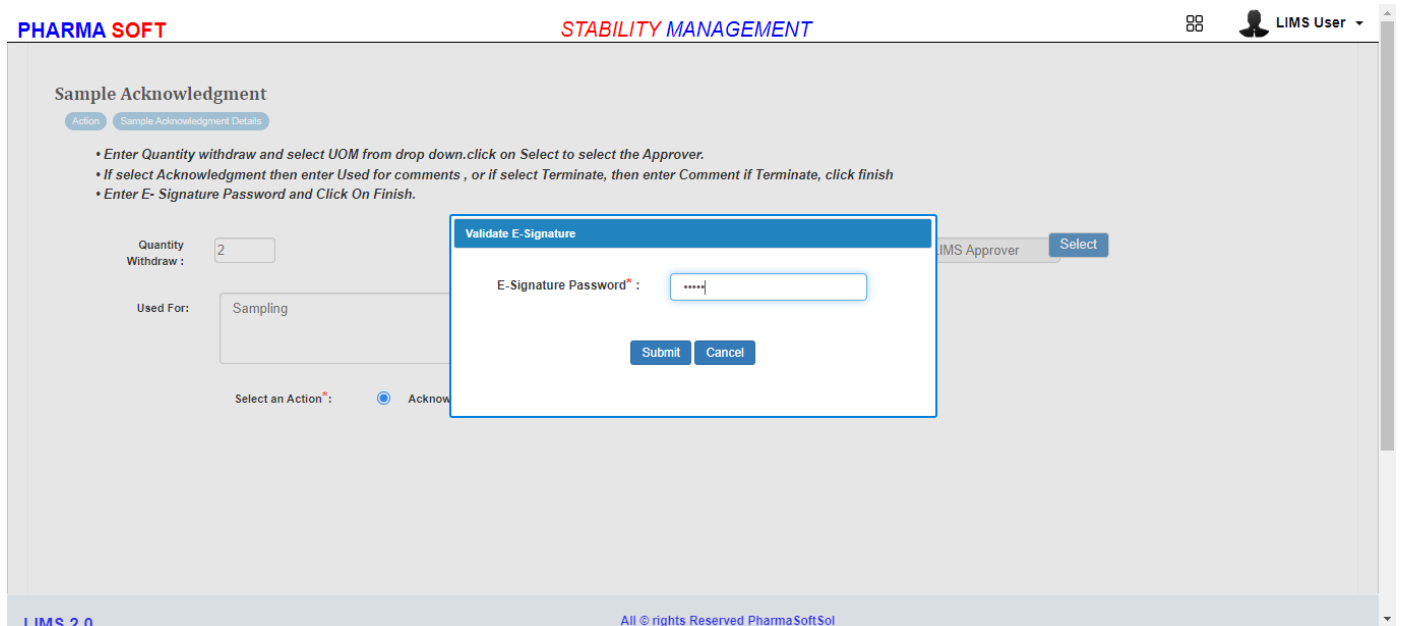


The screenshot shows the 'Sample Acknowledgment' form in the LIMS system. The form includes fields for 'Quantity Withdraw' (set to 2), 'Used For' (set to Sampling), and 'UOM' (set to Gm). A 'Validate E-Signature' modal is open, prompting for an 'E-Signature Password*'. The background form also shows an 'Approval From' dropdown set to 'LIMS Approver' and a 'Select' button. Instructions at the top of the form state: 'Enter Quantity withdraw and select UOM from drop down. click on Select to select the Approver.', 'If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish', and 'Enter E- Signature Password and Click On Finish.'.

12. Click on Submit

SSNo : PSS-LIMS-SS-012

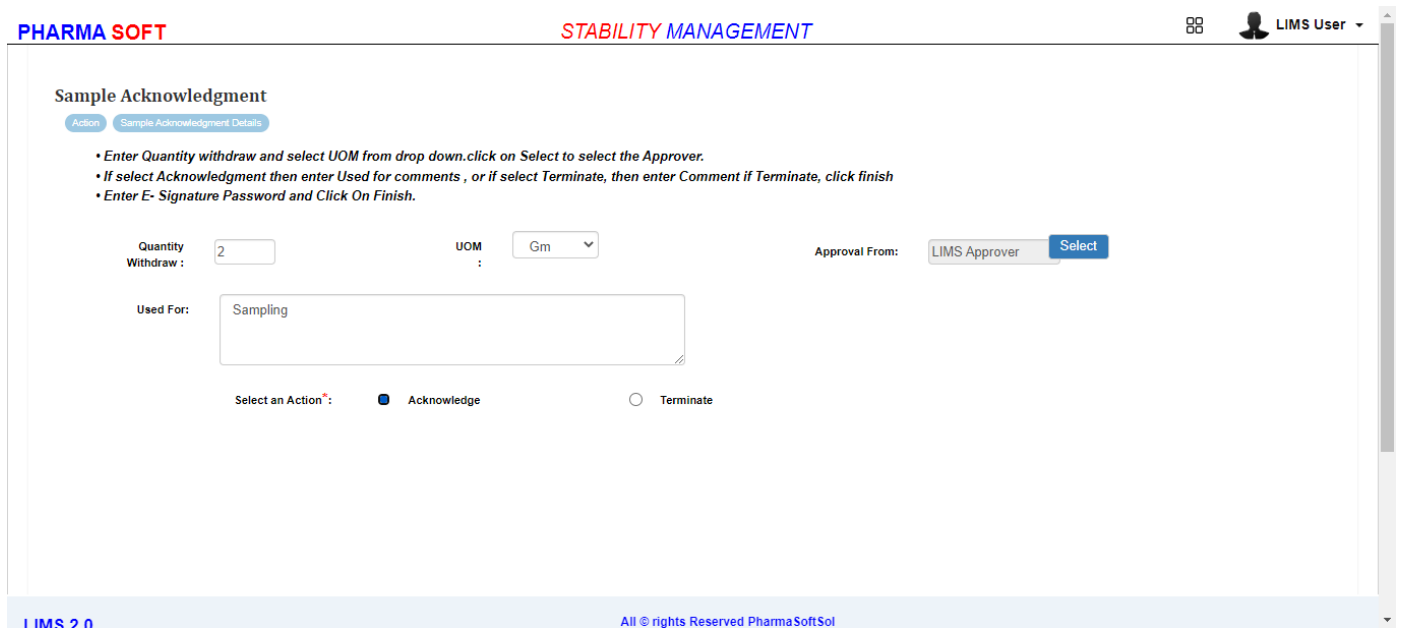

The screenshot shows the 'Sample Acknowledgment' form after the password validation step. The 'Quantity Withdraw' is 2, 'UOM' is Gm, and 'Used For' is Sampling. The 'Approval From' dropdown is set to 'LIMS Approver'. The 'Select an Action*' section shows the 'Acknowledge' radio button selected, with 'Terminate' as an alternative option. The 'Submit' button is visible at the bottom right of the form area. The same instructions from the previous step are present at the top.

	Test Case Name	ReSubmit RejectedSampleAcknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

13.Click on OK button

SSNo : PSS-LIMS-SS-013

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Sample Acknowledgment

Action **Sample Acknowledgment Details**

- Enter Quantity withdraw and select UOM from drop down.click on Select to select the Approver.
- If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish
- Enter E- Signature Password and Click On Finish.

Quantity Withdraw :
UOM :

Used For:

Select an Action*: ☒ Acknowledge

Success

Submitted successfully.

OK

LIMS Approver

LIMS 2.0
All © rights Reserved PharmaSoftSol

14.Click on User Name

SSNo : PSS-LIMS-SS-014

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Sample Acknowledgment

Action **Sample Acknowledgment Details**


- . Select Sample Acknowledgment, Click On Search, Select The Record and Click on Next.
- . If Reverted Sample Acknowledgments by Approver is Required, Then Select Reverted Sample Acknowledgments by Approver.

☐ Sample Acknowledgment
☒ Reverted Sample Acknowledgments by Approver

S.No	Product Name	Specification Number	AR Number	Protocol Number	Storage Conditions	Orientation	Stability Conditions	Packing Type	Batch Number	Start Date	RevertedBy/On	Comments
1	PARA65	SPSSPARA65	ASS220051	PC0007	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	Blister	Product Batch No 001	17/08/2022	LIMS User 17/08/2022	<input type="button" value="view"/>


Go to page: Row count:
Showing 1-1 of 1

LIMS 2.0
All © rights Reserved PharmaSoftSol


	Test Case Name	ReSubmit RejectedSampleAcknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

15.Click on Logout

SSNo : PSS-LIMS-SS-015

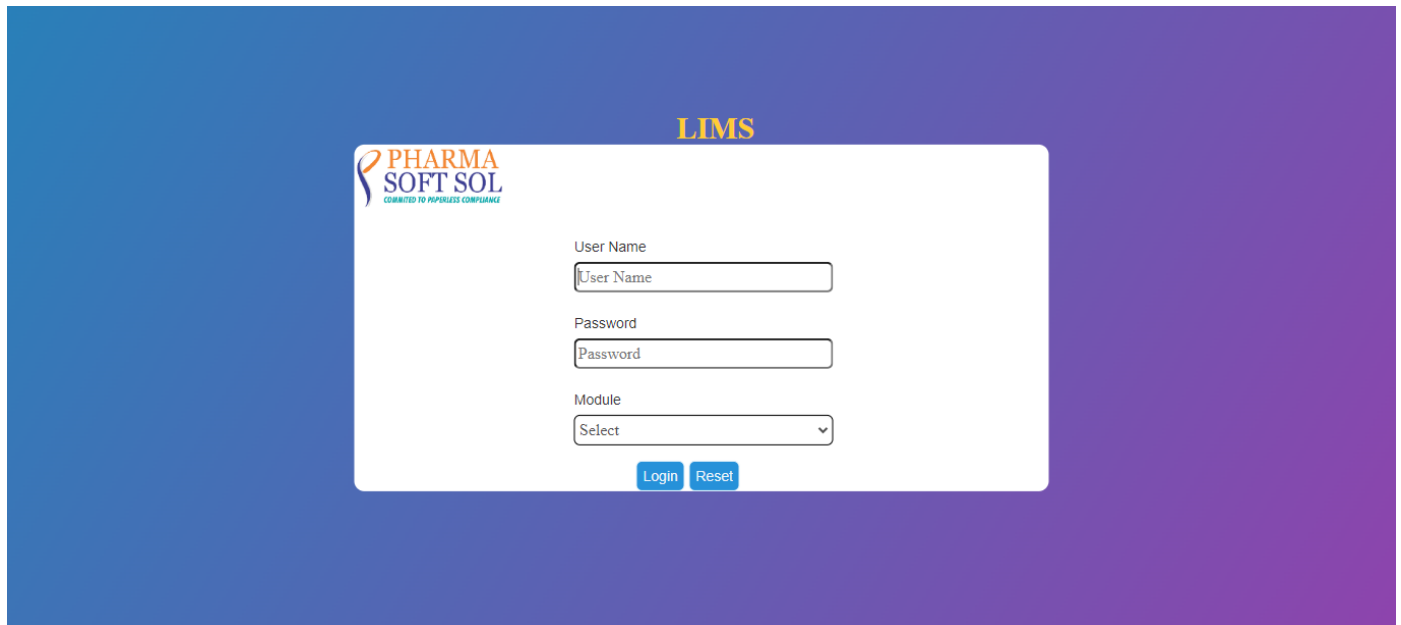

LIMS



User Name

Password

Module