	Test Case Name	Approve Sample Acknowledgement
	Test Case Number	LIMS-SM-014
	Actual Result	Pass

1. Enter the username, password, Select Module and click on login button

SSNo : PSS-LIMS-SS-001





2. Click on Sample Acknowledgement

SSNo : PSS-LIMS-SS-002

PHARMA SOFT

STABILITY MANAGEMENT



LIMS Approver

Sample Acknowledgment Approval

Activity :
 . Select The Record and Click On view.
 . Enter Comments and Click on Approve / Revert.
 . Enter E-Signature Password and Click on Submit.

S.No	Protocol Number	Product Name	Specification Number	AR Number	Storage Condition	Stability Condition	Packing type	Batch Number	Start Date	Done By/On	Action
1	PC0006	PARA65	SPSSPARA65	ASS220042	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	Blister	Product Batch No	16/08/2022	LIMS User /16/08/2022 19:15:55	View
2	PC0007	PARA65	SPSSPARA65	ASS220050	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	Blister	Product Batch No 001	17/08/2022	LIMS Approver /17/08/2022 14:53:43	View
3	PC0018	PARA65	SPPARACH001	ASS220081	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	Blister	Batch001	01/09/2022	LIMS User /01/09/2022 15:55:01	View

<< < 1 > >>


Go to page:

Row count:

Showing 1-3 of 3

LIMS 2.0

All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Approve Sample Acknowledgement
	Test Case Number	LIMS-SM-014
	Actual Result	Pass

3. Select a Record And Click on View

SSNo : PSS-LIMS-SS-003

PHARMA SOFT
STABILITY MANAGEMENT

 LIMS Approver

Sample Acknowledgement
Activity :
. Select The
. Enter Com
. Enter E-Sig
List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

Protocol Number:
Product Name:
Specification Number:
AR Number:
Batch Number:
Packing Type:
Storage Conditions:
Interval:
Qty Withdraw: UOM:
Total Qty Remained: UOM:
Start Date:
Used For:
Comments*:

View
View
View
Showing 1-3 of 3

LIMS 2.0

4. Enter Comments

SSNo : PSS-LIMS-SS-004

PHARMA SOFT
STABILITY MANAGEMENT

 LIMS Approver


Sample Acknowledgement
Activity :
. Select The
. Enter Com
. Enter E-Sig
List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

Protocol Number:
Product Name:
Specification Number:
AR Number:
Batch Number:
Packing Type:
Storage Conditions:
Interval:
Qty Withdraw: UOM:
Total Qty Remained: UOM:
Start Date:
Used For:
Comments*:

View
View
View
Showing 1-3 of 3

LIMS 2.0

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Approve Sample Acknowledgement
	Test Case Number	LIMS-SM-014
	Actual Result	Pass

5. Click on Approve

SSNo : PSS-LIMS-SS-005

PHARMA SOFT **STABILITY MANAGEMENT** ☰ 👤 LIMS Approver

Sample Acknowledgement

Activity :
 . Select The
 . Enter Com
 . Enter E-Sig

List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

<< < 1 > >>

Protocol Number: PC0018

Product Name: PARA65

Specification Number: SPPARACH001

AR Number: ASS220081

Batch Number: Batch001

Packing Type: Blister

Storage Conditions: 40°C ± 2°C / 7

Qty Withdraw: 2

Start Date: 01/09/2022

Used For: Sampling

Comments*: Approved

UOM: Gm

Approve Revert

Validate E-Signature

E-Signature Password* :

Submit Cancel

Action

View

View

View

Showing 1-3 of 3

LIMS 2.0

6. Enter E-Signature Password

SSNo : PSS-LIMS-SS-006

PHARMA SOFT **STABILITY MANAGEMENT** ☰ 👤 LIMS Approver

Sample Acknowledgement

Activity :
 . Select The
 . Enter Com
 . Enter E-Sig

List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

<< < 1 > >>

Protocol Number: PC0018

Product Name: PARA65

Specification Number: SPPARACH001

AR Number: ASS220081

Batch Number: Batch001

Packing Type: Blister

Storage Conditions: 40°C ± 2°C / 7

Qty Withdraw: 2

Start Date: 01/09/2022

Used For: Sampling

Comments*: Approved

UOM: Gm

Approve Revert

Validate E-Signature

E-Signature Password* :

Submit Cancel

Action


View

View

View


Showing 1-3 of 3

LIMS 2.0

	Test Case Name	Approve Sample Acknowledgement
	Test Case Number	LIMS-SM-014
	Actual Result	Pass

7. Click on Submit

SSNo : PSS-LIMS-SS-007

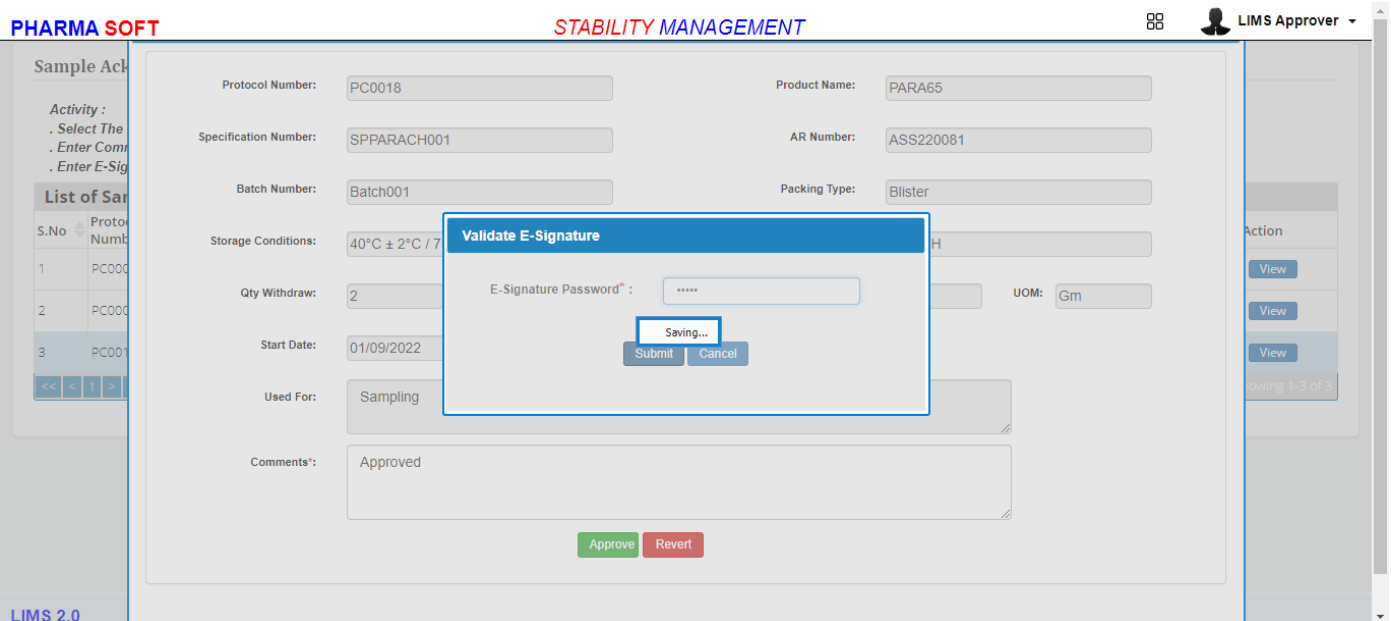

PHARMA SOFT **STABILITY MANAGEMENT** LIMS Approver

Sample Acknowledgement

Activity :
 . Select The
 . Enter Com
 . Enter E-Sig

List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

Protocol Number: PC0018 Product Name: PARA65

Specification Number: SPPARACH001 AR Number: ASS220081

Batch Number: Batch001 Packing Type: Blister

Storage Conditions: 40°C ± 2°C / 75%RH ± 5% MONTH

Qty Withdraw: 2 UOM: Gm

Start Date: 01/09/2022

Used For: Sampling

Comments*: Approved

Validate E-Signature

E-Signature Password*: *****

Saving...


Submit Cancel

Approve Revert

LIMS 2.0

8. Click on OK button

SSNo : PSS-LIMS-SS-008

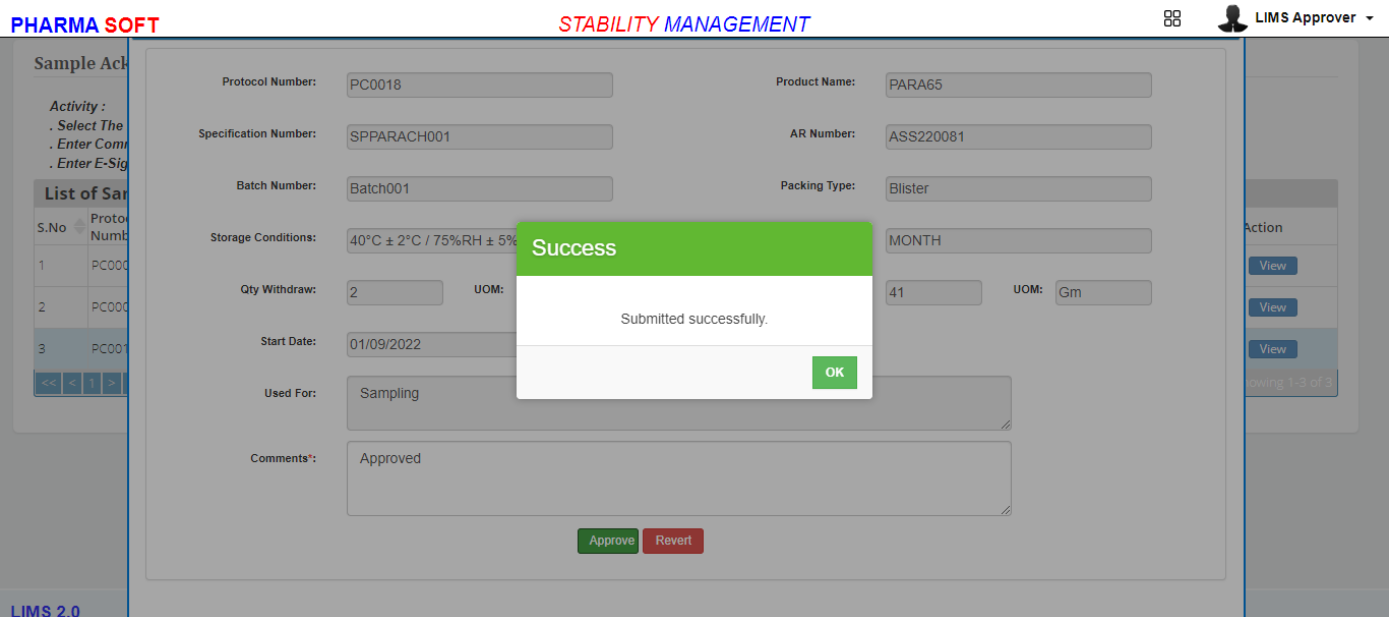

PHARMA SOFT **STABILITY MANAGEMENT** LIMS Approver

Sample Acknowledgement

Activity :
 . Select The
 . Enter Com
 . Enter E-Sig

List of Samples

S.No	Protocol Number
1	PC000
2	PC000
3	PC001

Protocol Number: PC0018 Product Name: PARA65

Specification Number: SPPARACH001 AR Number: ASS220081

Batch Number: Batch001 Packing Type: Blister

Storage Conditions: 40°C ± 2°C / 75%RH ± 5% MONTH

Qty Withdraw: 2 UOM: Gm

Start Date: 01/09/2022

Used For: Sampling

Comments*: Approved


Success

Submitted successfully.

OK

Approve Revert

LIMS 2.0


	Test Case Name	Approve Sample Acknowledgement
	Test Case Number	LIMS-SM-014
	Actual Result	Pass

9.Click on User Name

SSNo : PSS-LIMS-SS-009

PHARMA **SOFT SOL**
COMMITTED TO PAPERLESS COMPLIANCE

STABILITY MANAGEMENT

 LIMS Approver

PROTOCOL , MASTER OF CHAMBER & STORAGE CONDITION REGISTRATION

Create/Modify Protocol

Master Of Chamber

0 Protocol Approval

0 Master Of Chamber Approval

Recommendation

Status Logout

SCHEDULE

Schedule

Adhoc Schedule

0 Schedule Approval

0 Adhoc Schedule Approval

LOADING & SAMPLE ACKNOWLEDGMENT

Loading

Sample Acknowledgement


0 Loading Approval

2 Sample Acknowledgement Appr...

LIMS 2.0
All © rights Reserved PharmaSoftSol

10.Click on Logout

SSNo : PSS-LIMS-SS-010



LIMS

User Name

Password

Module

Select

Login

Reset