
	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

1.Enter the username, password, Select Module and click on login button


SSNo : PSS-LIMS-SS-001

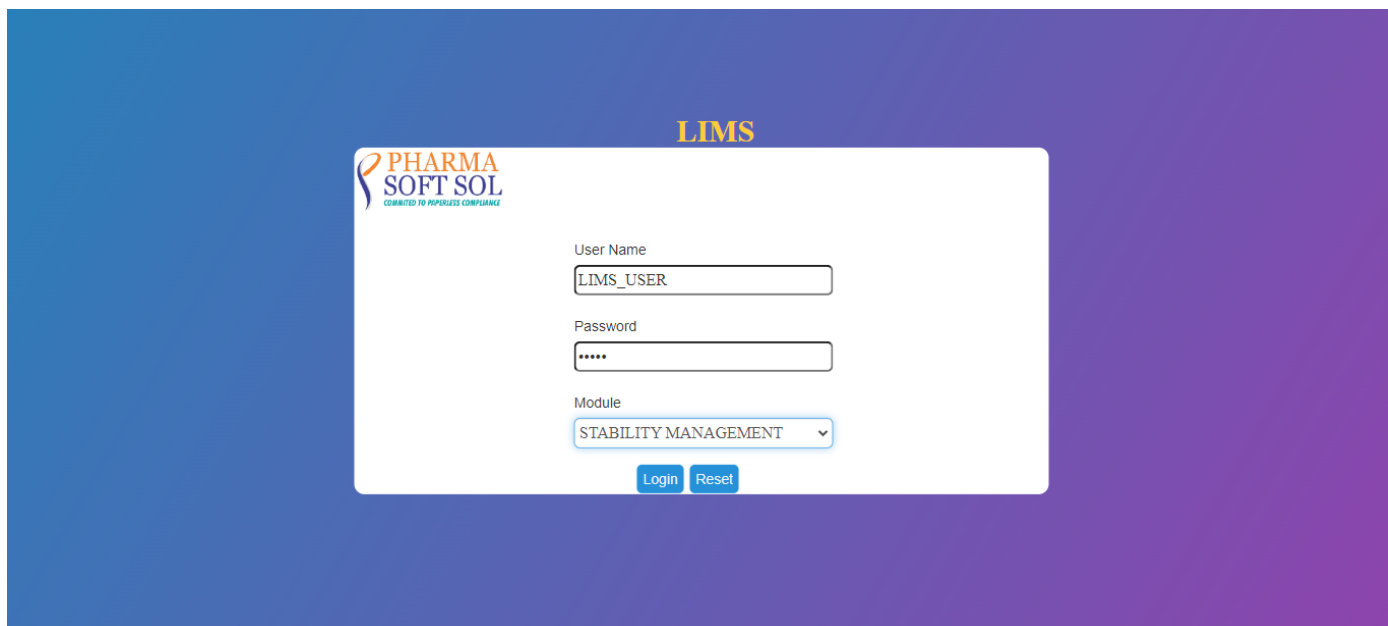

2.Click on Create

SSNo : PSS-LIMS-SS-002



	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

### 3. Click on Rejected

SSNo : PSS-LIMS-SS-003

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

**Create/Modify Protocol**

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

**Activity :**

- To create a protocol, select **Create Protocol** radio button and click on next.
- To Re-Submit the reverted protocols, select **Reverted Protocols by Approver** radio button. And select the required record in **List of Reverted Protocols by Approver** and click on next.
- To modify the protocol, select **Modify Existing Protocols** radio button and enter required search fields and click on search. And Select the required record in **List of Existed Protocols** and click on next.

☐ Create Protocol
☒ Reverted Protocols by Approver
☐ Modify Existing Protocols

☐ Save As Draft
☐ Issuance

**List of Reverted Protocols by Approver**

Protocol Number	Product Name	Generic Name	Specification Name	Reverted By/On	Comments
PC0014	PARA65	Paracetamol 650 mg	SPPARACHEMICAL	LIMS Approver 31/08/2022 18:21:40	<a href="#">View</a>

<<
<
1
>
>>
Go to page: 
Row count: 
Showing 1-1 of 1

IMS 2.0
All © rights Reserved PharmaSoftSol

### 4. Select a Record

SSNo : PSS-LIMS-SS-004

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

**Create/Modify Protocol**

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

**Activity :**

- To create a protocol, select **Create Protocol** radio button and click on next.
- To Re-Submit the reverted protocols, select **Reverted Protocols by Approver** radio button. And select the required record in **List of Reverted Protocols by Approver** and click on next.
- To modify the protocol, select **Modify Existing Protocols** radio button and enter required search fields and click on search. And Select the required record in **List of Existed Protocols** and click on next.

☐ Create Protocol
☒ Reverted Protocols by Approver
☐ Modify Existing Protocols


☐ Save As Draft
☐ Issuance

**List of Reverted Protocols by Approver**

Protocol Number	Product Name	Generic Name	Specification Name	Reverted By/On	Comments
PC0014	PARA65	Paracetamol 650 mg	SPPARACHEMICAL	LIMS Approver 31/08/2022 18:21:40	<a href="#">View</a>

<<
<
1
>
>>
Go to page: 
Row count: 
Showing 1-1 of 1

IMS 2.0
All © rights Reserved PharmaSoftSol

	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

5.Click on Next  
SSNo : PSS-LIMS-SS-005

PHARMA SOFT

STABILITY MANAGEMENT

LIMS User

Create/Modify Protocol

Action

Protocol Details

Batches & Product Information

Primary & Secondary Packing Material

Storage Conditions

Sample Details

Schedule Storage Condition

Modify Protocol Schedule

Schedule Modify

• Enter Purpose,Product Information,Product Description and Selection of Batches.

Objective\*

Formats

B

I

rotocol Objective

p

Responsibilities\*

Formats

B

I

rotocol Responsibilities

p

IMS 2.0

All © rights Reserved PharmaSoftSol

6.Clear Objective  
SSNo : PSS-LIMS-SS-006

PHARMA SOFT

STABILITY MANAGEMENT

LIMS User

Create/Modify Protocol

Action

Protocol Details

Batches & Product Information

Primary & Secondary Packing Material

Storage Conditions

Sample Details

Schedule Storage Condition

Modify Protocol Schedule

Schedule Modify

• Enter Purpose,Product Information,Product Description and Selection of Batches.

Objective\*

Formats

B

I

p

Responsibilities\*

Formats

B


I

rotocol Responsibilities

p

IMS 2.0

All © rights Reserved PharmaSoftSol

 <b>PHARMA SOFT SOL</b> <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

## 7.Enter Objective

SSNo : PSS-LIMS-SS-007

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

### Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter Purpose, Product Information, Product Description and Selection of Batches.

**Objective\***:

Formats
B
I
Text alignment icons
List icons

Protocol Objective

p

**Responsibilities\***:

Formats
B
I
Text alignment icons
List icons

rotocol Responsibilities

p

IMS 2.0
All © rights Reserved PharmaSoftSol

## 8.Clear Responsibilities

SSNo : PSS-LIMS-SS-008

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

### Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter Purpose, Product Information, Product Description and Selection of Batches.

**Objective\***:

Formats
B
I
Text alignment icons
List icons

Protocol Objective


p

**Responsibilities\***:

Formats
B
I
Text alignment icons
List icons

p

IMS 2.0
All © rights Reserved PharmaSoftSol

	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

9.Enter Responsibilities

SSNo : PSS-LIMS-SS-009

PHARMA SOFT

STABILITY MANAGEMENT

LIMS User

Create/Modify Protocol

Action

Protocol Details

Batches & Product Information

Primary & Secondary Packing Material

Storage Conditions

Sample Details

Schedule Storage Condition

Modify Protocol Schedule

Schedule Modify

• Enter Purpose,Product Information,Product Description and Selection of Batches.

Objective\*

Formats

B

I

Protocol Objective

p

Responsibilities\*

Formats

B

I

Protocol Responsibilities

p

IMS 2.0

All © rights Reserved PharmaSoftSol

10.Clear Stability Plan

SSNo : PSS-LIMS-SS-010

PHARMA SOFT

STABILITY MANAGEMENT

LIMS User

Create/Modify Protocol

Action

Protocol Details

Batches & Product Information

Primary & Secondary Packing Material

Storage Conditions

Sample Details

Schedule Storage Condition

Modify Protocol Schedule

Schedule Modify

• Enter Purpose,Product Information,Product Description and Selection of Batches.

Objective\*

Formats

B

I

Protocol Objective

p

Responsibilities\*

Formats

B

I


Protocol Responsibilities

p

IMS 2.0

All © rights Reserved PharmaSoftSol



	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

13.Enter Procedure

SSNo : PSS-LIMS-SS-013

PHARMA SOFT

STABILITY MANAGEMENT

LIMS User

p

Responsibilities\*

Formats

B

I

Protocol Responsibilities

p

Stability Plan\*

Formats

B

I

Protocol Stability Plan

p

Procedure\*

Formats

B

I

IMS 2.0

All © rights Reserved PharmaSoftSol

14.Click on Next

SSNo : PSS-LIMS-SS-014

PHARMA SOFT

STABILITY MANAGEMENT

LIMS User

Create/Modify Protocol

Action

Protocol Details

Batches & Product information

Primary & Secondary Packing Material

Storage Conditions

Sample Details


Schedule Storage Condition

Modify Protocol Schedule

Schedule Modify

IMS 2.0

All © rights Reserved PharmaSoftSol

 <b>PHARMA SOFT SOL</b> <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

## 15.Clear Manufacturing Date

SSNo : PSS-LIMS-SS-015

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter Number of Batches and click on Add,Enter Batch Details.

Number Of Batches\*:  Add

Batches					
S.No	Product Batch No	Mfg.Date	Exp.Date	Batch Size	Customer
1	44		08/31/2027	100	soft sol

Product Information:

• Enter Number of Rows for Production Information and click on Add.  
• Enter Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate Details and click on next.

No.of Rows for Product Information\*:  Add

Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate									
S.No	Product Batch No	API/PFI Batch no.	Name Of Material Used	AR Number	Mfg.Date	Exp.Date	Mfg.By	Is API	Is PFI
1	544	13	material	PSS-01	23-11-2020	23-12-2021	Pharma	<input checked="" type="checkbox"/>	<input type="checkbox"/>

IMS 2.0
All © rights Reserved PharmaSoftSol

## 16.Enter Manufacturing Date

SSNo : PSS-LIMS-SS-016

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter Number of Batches and click on Add,Enter Batch Details.

Number Of Batches\*:  Add

Batches					
S.No	Product Batch No	Mfg.Date	Exp.Date	Batch Size	Customer
1	44	08/31/2022	08/31/2027	100	soft sol

Product Information:


• Enter Number of Rows for Production Information and click on Add.  
• Enter Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate Details and click on next.

No.of Rows for Product Information\*:  Add

Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate									
S.No	Product Batch No	API/PFI Batch no.	Name Of Material Used	AR Number	Mfg.Date	Exp.Date	Mfg.By	Is API	Is PFI
1	544	13	material	PSS-01	23-11-2020	23-12-2021	Pharma	<input checked="" type="checkbox"/>	<input type="checkbox"/>

IMS 2.0
All © rights Reserved PharmaSoftSol



	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

## 17. Clear Expiry Date

SSNo : PSS-LIMS-SS-017

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter Number of Batches and click on Add, Enter Batch Details.

Number Of Batches\*:  Add

Batches					
S.No	Product Batch No	Mfg.Date	Exp.Date	Batch Size	Customer
1	44	08/31/2022		100	soft sol

Product Information:

• Enter Number of Rows for Production Information and click on Add.  
• Enter Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate Details and click on next.

No.of Rows for Product Information\*:  Add

Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate									
S.No	Product Batch No	API/PFI Batch no.	Name Of Material Used	AR Number	Mfg.Date	Exp.Date	Mfg.By	Is API	Is PFI
1	544	13	material	PSS-01	23-11-2020	23-12-2021	Pharma	<input checked="" type="checkbox"/>	<input type="checkbox"/>

IMS 2.0
All © rights Reserved PharmaSoftSol

## 18. Enter Expiry Date

SSNo : PSS-LIMS-SS-018

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter Number of Batches and click on Add, Enter Batch Details.

Number Of Batches\*:  Add

Batches					
S.No	Product Batch No	Mfg.Date	Exp.Date	Batch Size	Customer
1	44	08/31/2022	08/31/2028	100	soft sol


Product Information:

• Enter Number of Rows for Production Information and click on Add.  
• Enter Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate Details and click on next.

No.of Rows for Product Information\*:  Add

Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate									
S.No	Product Batch No	API/PFI Batch no.	Name Of Material Used	AR Number	Mfg.Date	Exp.Date	Mfg.By	Is API	Is PFI
1	544	13	material	PSS-01	23-11-2020	23-12-2021	Pharma	<input checked="" type="checkbox"/>	<input type="checkbox"/>

IMS 2.0
All © rights Reserved PharmaSoftSol

 <b>PHARMA SOFT SOL</b> <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

## 19. Clear Manufacturer Date

SSNo : PSS-LIMS-SS-019

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter Number of Batches and click on Add, Enter Batch Details.

Number Of Batches\*:  Add

Batches					
S.No	Product Batch No	Mfg.Date	Exp.Date	Batch Size	Customer
1	44	08/31/2022	08/31/2028	100	soft sol

Product Information:

• Enter Number of Rows for Production Information and click on Add.  
• Enter Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate Details and click on next.

No.of Rows for Product Information\*:  Add

Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate									
S.No	Product Batch No	API/PFI Batch no.	Name Of Material Used	AR Number	Mfg.Date	Exp.Date	Mfg.By	Is API	Is PFI
1	544	13	material	PSS-01		23-12-2021	Pharma	<input checked="" type="checkbox"/>	<input type="checkbox"/>

IMS 2.0
All © rights Reserved PharmaSoftSol

## 20. Enter Manufacturer Date

SSNo : PSS-LIMS-SS-020

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter Number of Batches and click on Add, Enter Batch Details.

Number Of Batches\*:  Add

Batches					
S.No	Product Batch No	Mfg.Date	Exp.Date	Batch Size	Customer
1	44	08/31/2022	08/31/2028	100	soft sol


Product Information:

• Enter Number of Rows for Production Information and click on Add.  
• Enter Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate Details and click on next.

No.of Rows for Product Information\*:  Add

Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate									
S.No	Product Batch No	API/PFI Batch no.	Name Of Material Used	AR Number	Mfg.Date	Exp.Date	Mfg.By	Is API	Is PFI
1	544	13	material	PSS-01	08/31/2022	23-12-2021	Pharma	<input checked="" type="checkbox"/>	<input type="checkbox"/>

IMS 2.0
All © rights Reserved PharmaSoftSol

 <b>PHARMA SOFT SOL</b> <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

## 21. Clear Expiry Date

SSNo : PSS-LIMS-SS-021

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter Number of Batches and click on Add, Enter Batch Details.

Number Of Batches\*:  Add

Batches					
S.No	Product Batch No	Mfg.Date	Exp.Date	Batch Size	Customer
1	44	08/31/2022	08/31/2028	100	soft sol

Product Information:

• Enter Number of Rows for Production Information and click on Add.  
• Enter Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate Details and click on next.

No.of Rows for Product Information\*:  Add

Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate									
S.No	Product Batch No	API/PFI Batch no.	Name Of Material Used	AR Number	Mfg.Date	Exp.Date	Mfg.By	Is API	Is PFI
1	544	13	material	PSS-01	08/31/2022		Pharma	<input checked="" type="checkbox"/>	<input type="checkbox"/>

IMS 2.0
All © rights Reserved PharmaSoftSol

## 22. Enter Expiry Date

SSNo : PSS-LIMS-SS-022

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter Number of Batches and click on Add, Enter Batch Details.

Number Of Batches\*:  Add

Batches					
S.No	Product Batch No	Mfg.Date	Exp.Date	Batch Size	Customer
1	44	08/31/2022	08/31/2028	100	soft sol


Product Information:

• Enter Number of Rows for Production Information and click on Add.  
• Enter Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate Details and click on next.

No.of Rows for Product Information\*:  Add

Active Pharmaceutical Ingredient/Pharmaceutical Formulation Intermediate									
S.No	Product Batch No	API/PFI Batch no.	Name Of Material Used	AR Number	Mfg.Date	Exp.Date	Mfg.By	Is API	Is PFI
1	544	13	material	PSS-01	08/31/2022	08/31/2034	Pharma	<input checked="" type="checkbox"/>	<input type="checkbox"/>

IMS 2.0
All © rights Reserved PharmaSoftSol

	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

23. Click on Next

SSNo : PSS-LIMS-SS-023

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

### Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter No. of Rows for Product Information and click on Add, Enter Primary Packing Material Details.

No. of Rows for Product Information\*:
Add

Primary Packing Material								
S.No	Product Batch No	Packing Component	Description	Type of Packing Material Used	Manufactured/Supplied By	Material Code	A.R.No/Ref.No	Qty per Container
1	10	round	2	Multiple types	biological	Resin	pp-01	3

• Enter No. of Rows for Product Information and click on Add, Enter Secondary Packing Material Details and click on next.

No. of Rows for Product Information:
Add

Secondary Packing Material								
S.No	Product Batch No	Packing Material	Description	Type of Packing Material Used	Manufactured/Supplied By	Material Code	A.R.No/Ref.No	Qty per Container

IMS 2.0
All © rights Reserved PharmaSoftSol

24. Click on Next

SSNo : PSS-LIMS-SS-024

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

### Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

Enter Primary Packing Material Details.


Type of Packing Material Used	Manufactured/Supplied By	Material Code	A.R.No/Ref.No	Qty per Container
Multiple types	biological	Resin	pp-01	3

Enter Secondary Packing Material Details and click on next.

Type of Packing Material Used	Manufactured/Supplied By	Material Code	A.R.No/Ref.No	Qty per Container
-------------------------------	--------------------------	---------------	---------------	-------------------

• click on Select to select the required Storage Conditions.

IMS 2.0
All © rights Reserved PharmaSoftSol

	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

## 25. Clear Sampling & Labeling

SSNo : PSS-LIMS-SS-025


**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• click on Select to select the required Storage Conditions. Enter Sampling & Labeling details and click on next.

Storage Conditions\*: Select

Proposed Condition of Stability Study List			
Storage Conditions	Sample Orientation	Stability Conditions	
40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	

Sampling & Labeling\*: 

Maximum allowed characters should be 1500

IMS 2.0

All © rights Reserved PharmaSoftSol

## 26. Enter Sampling & Labeling

SSNo : PSS-LIMS-SS-026


**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• click on Select to select the required Storage Conditions. Enter Sampling & Labeling details and click on next.

Storage Conditions\*: Select


Proposed Condition of Stability Study List			
Storage Conditions	Sample Orientation	Stability Conditions	
40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	

Sampling & Labeling\*: 

Protocol Created

IMS 2.0

All © rights Reserved PharmaSoftSol

	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

27. Click on Next


SSNo : PSS-LIMS-SS-027

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

Click on the Stability Conditions tab and click on next.

Stability Conditions	
ACCELERATED	

• Enter No. of Rows for Sample Details and click on Add.

IMS 2.0

All © rights Reserved PharmaSoftSol

28. Enter No. of Rows

SSNo : PSS-LIMS-SS-028

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter No. of Rows for Sample Details and click on Add.

No. of Rows For Sample Details\*:  Add

• Enter Sample Details Kept For Stability Details and click on Browse and select the required document. Enter Remarks and click on next.

Sample Details Kept For Stability																													
S.No	Storage Conditions	Orientation	Stability Conditions	Add Batch Details	Product Batch No	Mfg. Date	Exp. Date	Batch Size	Customer	Quantity Charged	UOM																		
<div> <b>Existed Sample Details Kept For Stability</b> <table border="1"> <thead> <tr> <th>Storage Condition</th> <th>Stability Condition</th> <th>Batch Number</th> <th>Batch Size</th> <th>Mfg Date</th> <th>Exp Date</th> <th>Customer</th> <th>Quantity Charged</th> <th>UOM</th> </tr> </thead> <tbody> <tr> <td>40°C ± 2°C / 75%RH ± 5%RH</td> <td>ACCELERATED</td> <td>44</td> <td>100</td> <td>08/31/2022</td> <td>08/31/2027</td> <td>soft sol</td> <td>3</td> <td>Gm</td> </tr> </tbody> </table> </div>												Storage Condition	Stability Condition	Batch Number	Batch Size	Mfg Date	Exp Date	Customer	Quantity Charged	UOM	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	44	100	08/31/2022	08/31/2027	soft sol	3	Gm
Storage Condition	Stability Condition	Batch Number	Batch Size	Mfg Date	Exp Date	Customer	Quantity Charged	UOM																					
40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	44	100	08/31/2022	08/31/2027	soft sol	3	Gm																					

<<
<
1
>
>>
Go to page:  Row count:  Showing 1-1 of 1

Existing Document:
View


Upload Document:
Browse

Remarks\*:

Protocol Created

IMS 2.0

All © rights Reserved PharmaSoftSol

	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

29. Click on Add

SSNo : PSS-LIMS-SS-029

**PHARMA SOFT**
**STABILITY MANAGEMENT**

LIMS User

Create/Modify Protocol

Action

Protocol Details

Batches & Product Information

Primary & Secondary Packing Material

Storage Conditions

Sample Details

Schedule Storage Condition

Modify Protocol Schedule

Schedule Modify

• Enter No. of Rows for Sample Details and click on Add.

No. of Rows For Sample Details\*:  Add

• Enter Sample Details Kept For Stability Details and click on Browse and select the required document. Enter Remarks and click on next.

Sample Details Kept For Stability

S.No	Storage Conditions	Orientation	Stability Conditions	Add Batch Details	Product Batch No	Mfg. Date	Exp. Date	Batch Size	Customer	Quantity Charged	UOM
1				<span>Add</span>							Se

Existed Sample Details Kept For Stability

Storage Condition	Stability Condition	Batch Number	Batch Size	Mfg Date	Exp Date	Customer	Quantity Charged	UOM
40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	44	100	08/31/2022	08/31/2027	soft sol	3	Gm

<< < 1 > >> Go to page:  Row count:  Showing 1-1 of 1

Existing Document:

SamplePDF.pdf

View

Upload Document:

Browse

Remarks\*

Protocol Created

IIMS 2.0

All © rights Reserved PharmaSoftSol

30. Click on Add

SSNo : PSS-LIMS-SS-030

**PHARMA SOFT**
**STABILITY MANAGEMENT**

LIMS User

Batch and Storage Details

Batch and Storage Condition Details

sn#	Product Batch No	Batch Size	Manufacturing Date	Expiry Date	Customer	Storage Condition	Orientation Name	Stability Condition
1	44	100	08/31/2022	08/31/2028	soft sol	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED


Select

IIMS 2.0

All © rights Reserved PharmaSoftSol

© SmartPharmaSoftSol

Page 15 of 21

 <b>PHARMA SOFT SOL</b> <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

31.select Batch

SSNo : PSS-LIMS-SS-031

**PHARMA SOFT** **STABILITY MANAGEMENT** ☰ LIMS User

Create/Modify Protocol

**Batch and Storage Details**

sno	Product Batch No	Batch Size	Manufacturing Date	Expiry Date	Customer	Storage Condition	Orientation Name	Stability Condition
1	44	100	08/31/2022	08/31/2028	soft sol	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED

[Select](#)

IMS 2.0 All © rights Reserved PharmaSoftSol

32.Click on Select Button

SSNo : PSS-LIMS-SS-032

**PHARMA SOFT** **STABILITY MANAGEMENT** ☰ LIMS User

Create/Modify Protocol

[Action](#) [Protocol Details](#) [Batches & Product Information](#) [Primary & Secondary Packing Material](#) [Storage Conditions](#) [Sample Details](#) [Schedule Storage Condition](#) [Modify Protocol Schedule](#) [Schedule Modify](#)

• Enter No.of Rows for Sample Details and click on Add.

No.of Rows For Sample Details\*:  [Add](#)

• Enter Sample Details Kept For Stability Details and click on Browse and select the required document.Enter Remarks and click on next.

**Sample Details Kept For Stability**

S.No	Storage Conditions	Orientation	Stability Conditions	Add Batch Details	Product Batch No	Mfg.Date	Exp.Date	Batch Size	Customer	Quantity Charged	UOM
1	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	<a href="#">Add</a>	44	08/31/2022	08/31/2028	100	soft sol	<input type="text"/>	<input type="text"/>

**Existed Sample Details Kept For Stability**


Storage Condition	Stability Condition	Batch Number	Batch Size	Mfg Date	Exp Date	Customer	Quantity Charged	UOM
40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	44	100	08/31/2022	08/31/2027	soft sol	3	Gm

<< < 1 > >> Go to page:  Row count:  Showing 1-1 of 1

Existing Document:  [View](#) Upload Document:  [Browse](#) Remarks\*:

IMS 2.0 All © rights Reserved PharmaSoftSol



	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

### 33. Enter Quantity Charged

SSNo : PSS-LIMS-SS-033

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter No. of Rows for Sample Details and click on Add.

No. of Rows For Sample Details\*:  Add

• Enter Sample Details Kept For Stability Details and click on Browse and select the required document. Enter Remarks and click on next.

S.No	Storage Conditions	Orientation	Stability Conditions	Add Batch Details	Product Batch No	Mfg. Date	Exp. Date	Batch Size	Customer	Quantity Charged	UOM	
1	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	<button>Add</button>	44	08/31/2022	08/31/2028	100	soft sol	<input type="text" value="3"/>	<input type="text" value="Gm"/>	

Storage Condition	Stability Condition	Batch Number	Batch Size	Mfg Date	Exp Date	Customer	Quantity Charged	UOM
40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	44	100	08/31/2022	08/31/2027	soft sol	3	Gm

<< < 1 > >> Go to page:  Row count:  Showing 1-1 of 1

Existing Document:  View Upload Document:  Browse Remarks\*:

IMS 2.0
All © rights Reserved PharmaSoftSol

### 34. Select UOM

SSNo : PSS-LIMS-SS-034

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action
Protocol Details
Batches & Product Information
Primary & Secondary Packing Material
Storage Conditions
Sample Details
Schedule Storage Condition
Modify Protocol Schedule
Schedule Modify

• Enter No. of Rows for Sample Details and click on Add.

No. of Rows For Sample Details\*:  Add

• Enter Sample Details Kept For Stability Details and click on Browse and select the required document. Enter Remarks and click on next.


S.No	Storage Conditions	Orientation	Stability Conditions	Add Batch Details	Product Batch No	Mfg. Date	Exp. Date	Batch Size	Customer	Quantity Charged	UOM	
1	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	<button>Add</button>	44	08/31/2022	08/31/2028	100	soft sol	<input type="text" value="3"/>	<input type="text" value="G"/>	

Storage Condition	Stability Condition	Batch Number	Batch Size	Mfg Date	Exp Date	Customer	Quantity Charged	UOM
40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	44	100	08/31/2022	08/31/2027	soft sol	3	Gm

<< < 1 > >> Go to page:  Row count:  Showing 1-1 of 1

Existing Document:  View Upload Document:  Browse Remarks\*:

IMS 2.0
All © rights Reserved PharmaSoftSol

	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

### 35.Upload Document1

SSNo : PSS-LIMS-SS-035

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action Protocol Details Batches & Product Information Primary & Secondary Packing Material Storage Conditions **Sample Details** Schedule Storage Condition Modify Protocol Schedule Schedule Modify

• Enter No.of Rows for Sample Details and click on Add.

No.of Rows For Sample Details\*:  **Add**

• Enter Sample Details Kept For Stability Details and click on Browse and select the required document.Enter Remarks and click on next.

S.No	Storage Conditions	Orientation	Stability Conditions	Add Batch Details	Product Batch No	Mfg.Date	Exp.Date	Batch Size	Customer	Quantity Charged	UOM	
1	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	<b>Add</b>	44	08/31/2022	08/31/2028	100	soft sol	<input type="text" value="3"/>	<b>G</b>	

Storage Condition	Stability Condition	Batch Number	Batch Size	Mfg Date	Exp Date	Customer	Quantity Charged	UOM
40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	44	100	08/31/2022	08/31/2027	soft sol	3	Gm

<< < 1 > >> Go to page:  Row count:  Showing 1-1 of 1

Existing Document:  **View** Upload Document:  **Browse** Remarks\*:

IMS 2.0
All © rights Reserved PharmaSoftSol

### 36.Click on Next

SSNo : PSS-LIMS-SS-036

**PHARMA SOFT**
**STABILITY MANAGEMENT**
LIMS User

Create/Modify Protocol

Action Protocol Details Batches & Product Information Primary & Secondary Packing Material Storage Conditions **Sample Details** **Schedule Storage Condition** Modify Protocol Schedule Schedule Modify

• Enter Remarks and click on next.


Exp.Date	Batch Size	Customer	Quantity Charged	UOM	
08/31/2028	100	soft sol	<input type="text" value="3"/>	<b>G</b>	

Date	Customer	Quantity Charged	UOM
31/2027	soft sol	3	Gm

Showing 1-1 of 1


**se** Remarks\*:

IMS 2.0
All © rights Reserved PharmaSoftSol

	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

37. Click on Finish

SSNo : PSS-LIMS-SS-037

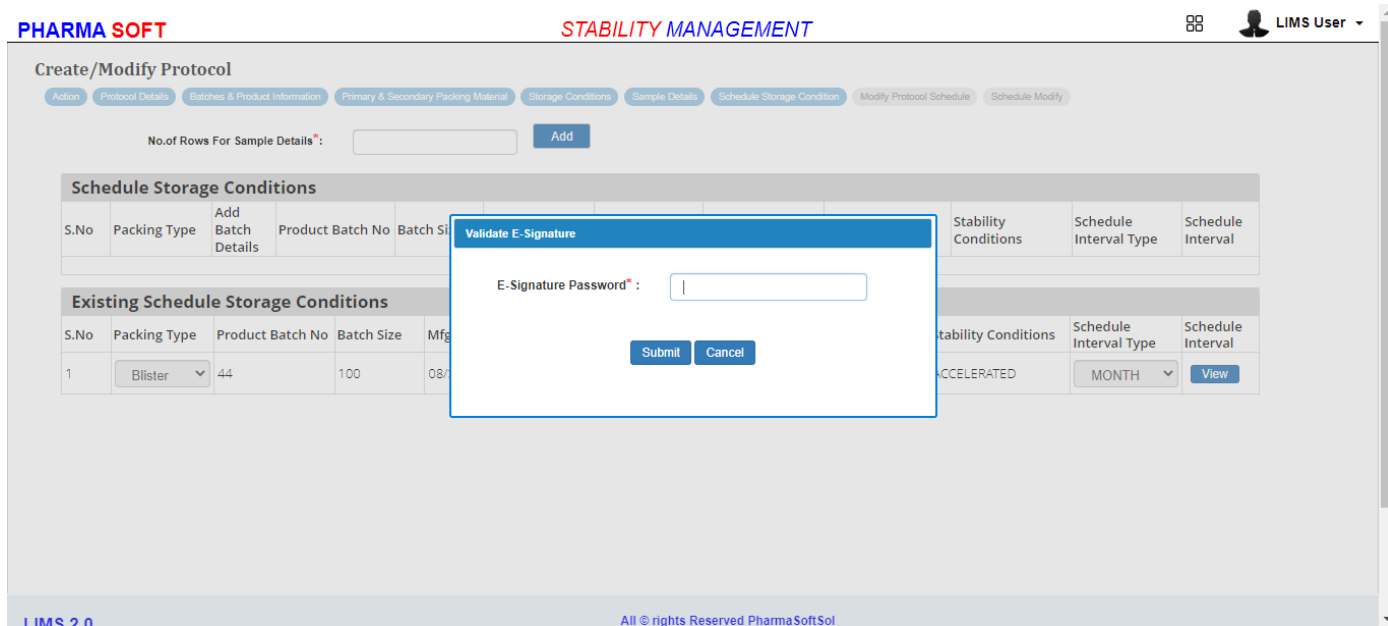

**PHARMA SOFT** **STABILITY MANAGEMENT** LIMS User

Create/Modify Protocol

Action Protocol Details Batches & Product Information Primary & Secondary Packing Material Storage Conditions Sample Details Schedule Storage Condition Modify Protocol Schedule Schedule Modify

No. of Rows For Sample Details\*:  Add

S.No	Packing Type	Add Batch Details	Product Batch No	Batch Size	Mfg	Stability Conditions	Schedule Interval Type	Schedule Interval
1	Blister	44	100	08/		ACCELERATED	MONTH	View

**Existing Schedule Storage Conditions**

S.No	Packing Type	Product Batch No	Batch Size	Mfg	Stability Conditions	Schedule Interval Type	Schedule Interval
1	Blister	44	100	08/	ACCELERATED	MONTH	View

Validate E-Signature


E-Signature Password\*:

Submit Cancel

LIMS 2.0 All © rights Reserved PharmaSoftSol

38. Enter E-Signature Password

SSNo : PSS-LIMS-SS-038

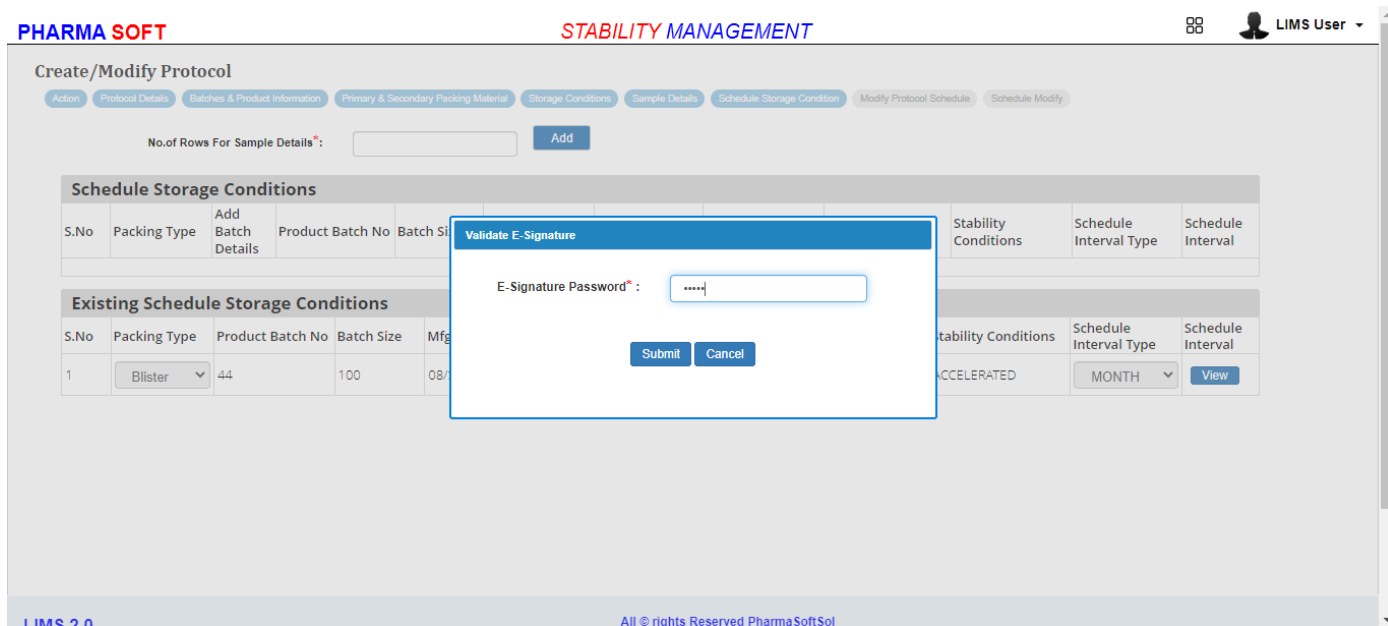

**PHARMA SOFT** **STABILITY MANAGEMENT** LIMS User

Create/Modify Protocol

Action Protocol Details Batches & Product Information Primary & Secondary Packing Material Storage Conditions Sample Details Schedule Storage Condition Modify Protocol Schedule Schedule Modify

No. of Rows For Sample Details\*:  Add

S.No	Packing Type	Add Batch Details	Product Batch No	Batch Size	Mfg	Stability Conditions	Schedule Interval Type	Schedule Interval
1	Blister	44	100	08/		ACCELERATED	MONTH	View

**Existing Schedule Storage Conditions**


S.No	Packing Type	Product Batch No	Batch Size	Mfg	Stability Conditions	Schedule Interval Type	Schedule Interval
1	Blister	44	100	08/	ACCELERATED	MONTH	View

Validate E-Signature

E-Signature Password\*:


Submit Cancel

LIMS 2.0 All © rights Reserved PharmaSoftSol

	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

39. Click on Submit

SSNo : PSS-LIMS-SS-039

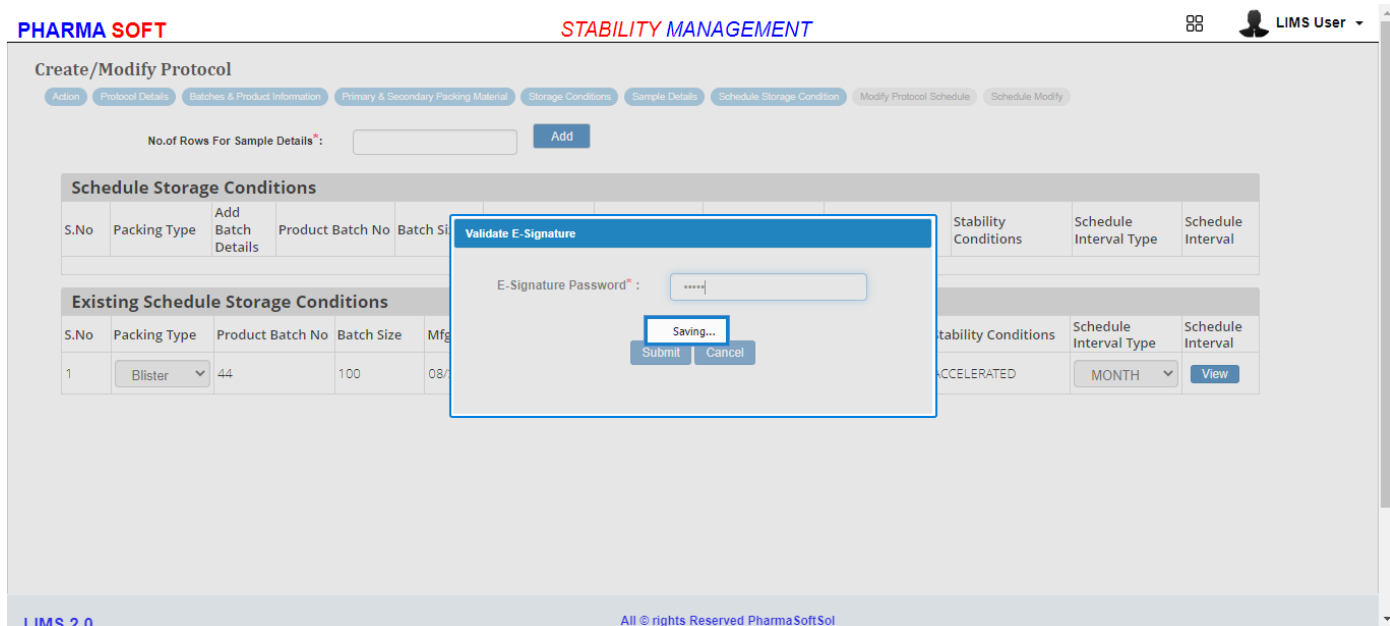

**PHARMA SOFT** **STABILITY MANAGEMENT** LIMS User

Create/Modify Protocol

[Action](#)
[Protocol Details](#)
[Batches & Product Information](#)
[Primary & Secondary Packing Material](#)
[Storage Conditions](#)
[Sample Details](#)
[Schedule Storage Condition](#)
[Modify Protocol Schedule](#)
[Schedule Modify](#)

No. of Rows For Sample Details\*:  [Add](#)

S.No	Packing Type	Add Batch Details	Product Batch No	Batch Size	Mfg. Date	Stability Conditions	Schedule Interval Type	Schedule Interval
1	Blister	44	100	08/31/2022	ACCELERATED	MONTH	<a href="#">View</a>	

**Existing Schedule Storage Conditions**

S.No	Packing Type	Product Batch No	Batch Size	Mfg. Date	Stability Conditions	Schedule Interval Type	Schedule Interval
1	Blister	44	100	08/31/2022	ACCELERATED	MONTH	<a href="#">View</a>

**Validate E-Signature**

E-Signature Password\*:


[Saving...](#)

[Submit](#) [Cancel](#)

LIMS 2.0 All © rights Reserved PharmaSoftSol

40. Click on OK button

SSNo : PSS-LIMS-SS-040

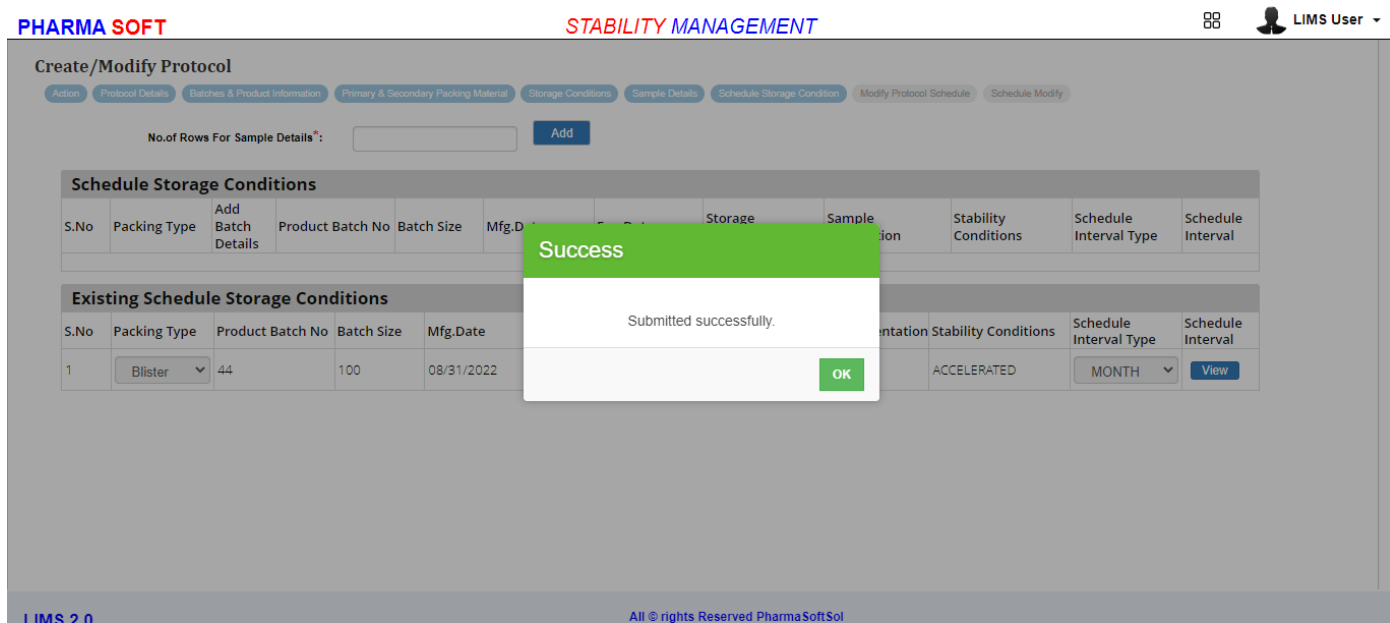

**PHARMA SOFT** **STABILITY MANAGEMENT** LIMS User

Create/Modify Protocol

[Action](#)
[Protocol Details](#)
[Batches & Product Information](#)
[Primary & Secondary Packing Material](#)
[Storage Conditions](#)
[Sample Details](#)
[Schedule Storage Condition](#)
[Modify Protocol Schedule](#)
[Schedule Modify](#)

No. of Rows For Sample Details\*:  [Add](#)

S.No	Packing Type	Add Batch Details	Product Batch No	Batch Size	Mfg. Date	Stability Conditions	Schedule Interval Type	Schedule Interval
1	Blister	44	100	08/31/2022	ACCELERATED	MONTH	<a href="#">View</a>	

**Existing Schedule Storage Conditions**


S.No	Packing Type	Product Batch No	Batch Size	Mfg. Date	Stability Conditions	Schedule Interval Type	Schedule Interval
1	Blister	44	100	08/31/2022	ACCELERATED	MONTH	<a href="#">View</a>

**Success**

Submitted successfully.

[OK](#)

LIMS 2.0 All © rights Reserved PharmaSoftSol

	<b>Test Case Name</b>	Reinitiate Rejected Protocol
	<b>Test Case Number</b>	PSS-LIMS-006
	<b>Actual Result</b>	Pass

41.Click on User Name

SSNo : PSS-LIMS-SS-041



42.Click on Logout

SSNo : PSS-LIMS-SS-042

