
	Test Case Name	Create Sample Acknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

1.Enter the username, password, Select Module and click on login button

SSNo : PSS-LIMS-SS-001

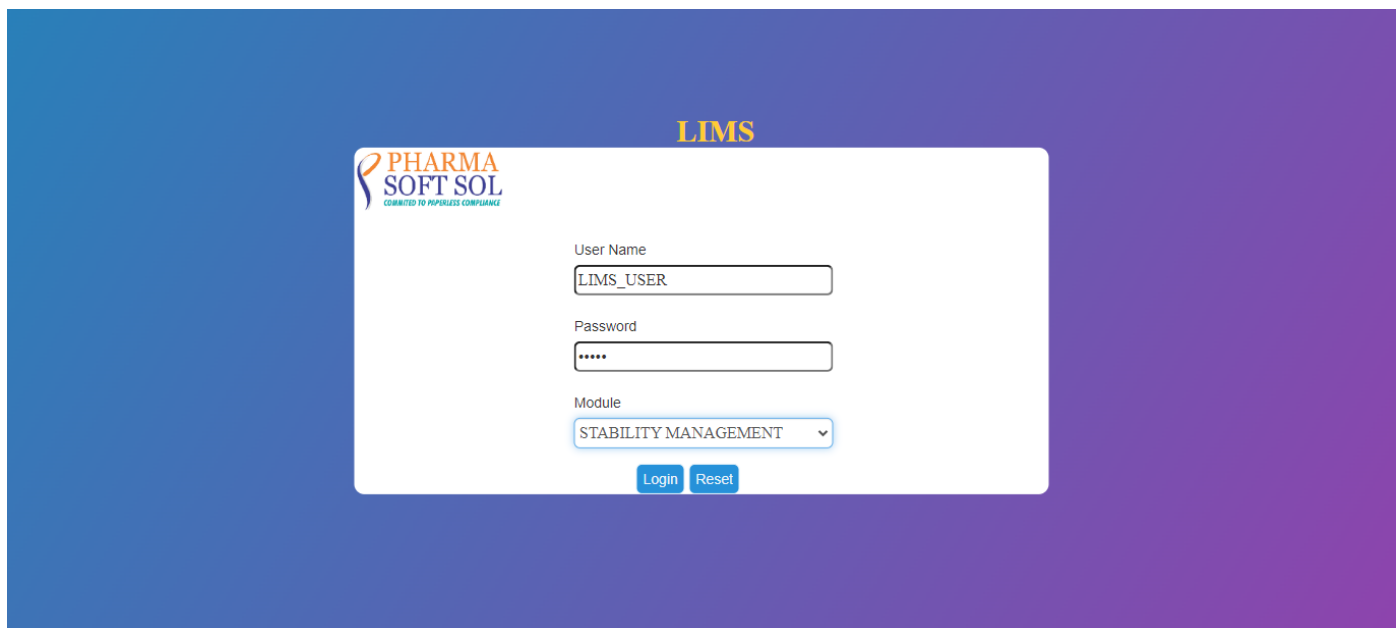

The screenshot shows the LIMS login interface. At the top, the 'LIMS' title is displayed in yellow. Below it, the PharmaSoftSol logo is visible. The login form contains three input fields: 'User Name' with the value 'LIMS_USER', 'Password' with masked characters '*****', and 'Module' with a dropdown menu set to 'STABILITY MANAGEMENT'. At the bottom of the form are two buttons: 'Login' and 'Reset'.

2.Click on Sample Acknowledgement


SSNo : PSS-LIMS-SS-002



The screenshot displays the 'Sample Acknowledgement' screen within the 'STABILITY MANAGEMENT' module. The header shows 'PHARMA SOFT' and 'STABILITY MANAGEMENT'. The user is logged in as 'LIMS USER'. The main content area has a tabbed interface with 'Sample Acknowledgement Details' selected. Below the tabs, there are instructions: '. Select Sample Acknowledgement, Click On Search, Select The Record and Click on Next.' and '. If Reverted Sample Acknowledgments by Approver is Required, Then Select Reverted Sample Acknowledgments by Approver.' Two radio buttons are present: 'Sample Acknowledgement' (unselected) and 'Reverted Sample Acknowledgments by Approver' (selected). Below the radio buttons is a table titled 'List of Reverted Sample Acknowledgements by Approver'.

S.No	Protocol Id	Product Name	Specification Number	AR Number	Protocol Number	Storage Conditions	Stability Conditions	Packing Type	Batch Number	Start Date	RevertedBy/On	Comments
1	7	Aspirin	SNSS 001	ASS220013	PS0003	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	Box	Product Batch No 1	11/06/2022	LIMS APPROVER 13/07/2022	view

Below the table, there is a pagination bar with navigation arrows, 'Go to page: 1', 'Row count: 10', and 'Showing 1-1 of 1'.

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Sample Acknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

3. Click on Sample Acknowledgment

SSNo : PSS-LIMS-SS-003

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Sample Acknowledgment

Action
Sample Acknowledgment Details

. Select Sample Acknowledgment, Click On Search, Select The Record and Click on Next.
 . If Reverted Sample Acknowledgments by Approver is Required, Then Select Reverted Sample Acknowledgments by Approver.

☒ Sample Acknowledgment
 ☐ Reverted Sample Acknowledgments by Approver

From Date :
 To Date:

Product : Browse
 Specification:

Search Clear

Samples Loaded

S.No	Product	Specification Number	AR Number	Protocol Number	Storage Condition	Stability Condition	Orientation	Packing Type	Batch Number	Start Date	Total Qty Remained	Qty Remained UOM	Interval	Is deviated
No data available!														

Go to page: Row count:

IMS 2.0
All © rights Reserved PharmaSoftSol

4. Click on Search

SSNo : PSS-LIMS-SS-004

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Sample Acknowledgment

Action
Sample Acknowledgment Details

. Select Sample Acknowledgment, Click On Search, Select The Record and Click on Next.
 . If Reverted Sample Acknowledgments by Approver is Required, Then Select Reverted Sample Acknowledgments by Approver.

☒ Sample Acknowledgment
 ☐ Reverted Sample Acknowledgments by Approver

From Date :
 To Date:

Product : Browse
 Specification:


Search Clear

Samples Loaded

S.No	Product	Specification Number	AR Number	Protocol Number	Storage Condition	Stability Condition	Orientation	Packing Type	Batch Number	Start Date	Total Qty Remained	Qty Remained UOM	Interval	Is deviated
No data available!														

Go to page: Row count:

IMS 2.0
All © rights Reserved PharmaSoftSol

	Test Case Name	Create Sample Acknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

5. Select a Record

SSNo : PSS-LIMS-SS-005

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Sample Acknowledgment

Action Sample Acknowledgment Details

. Select Sample Acknowledgment, Click On Search, Select The Record and Click on Next.
. If Reverted Sample Acknowledgments by Approver is Required, Then Select Reverted Sample Acknowledgments by Approver.

☒ Sample Acknowledgment
☐ Reverted Sample Acknowledgments by Approver

From Date: To Date:
Product: **Browse** Specification:

Search **Clear**

Samples Loaded														
S.No	Product	Specification Number	AR Number	Protocol Number	Storage Condition	Stability Condition	Orientation	Packing Type	Batch Number	Start Date	Total Qty. Remained	Qty Remained UOM	Interval	Is deviated
1	Aspirin	SNSS 001	ASS220014	PS0003	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	UP RIGHT	Box	Product Batch No 1	11/06/2022	43	Mg	3 MONTH	No
					40°C ± 2°C /				Product					

IMS 2.0
All © rights Reserved PharmaSoftSol

6. Click on Next

SSNo : PSS-LIMS-SS-006

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Sample Acknowledgment

Action Sample Acknowledgment Details

On Search, Select The Record and Click on Next.
y Approver is Required, Then Select Reverted Sample Acknowledgments by Approver.


☐ Reverted Sample Acknowledgments by Approver

To Date:
 Browse Specification:

Search **Clear**

Protocol Number	Storage Condition	Stability Condition	Orientation	Packing Type	Batch Number	Start Date	Total Qty. Remained	Qty Remained UOM	Interval	Is deviated
PS0003	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	UP RIGHT	Box	Product Batch No 1	11/06/2022	43	Mg	3 MONTH	No
	40°C ± 2°C /				Product					

IMS 2.0
All © rights Reserved PharmaSoftSol

	Test Case Name	Create Sample Acknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

7.Enter Quantity
SSNo : PSS-LIMS-SS-007

PHARMA SOFT

STABILITY MANAGEMENT

LIMS USER

Sample Acknowledgment

Action

Sample Acknowledgment Details

• Enter Quantity withdraw and select UOM from drop down.click on Select to select the Approver.

• If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish

• Enter E- Signature Password and Click On Finish.

Quantity Withdraw :

2

UOM :

Select

Approval From:

Select

Used For:

Maximum allowed characters should be 100

Select an Action*:

☐ Acknowledge

☐ Terminate

I IMS 2.0

All © rights Reserved PharmaSoftSol

8.Select UOM
SSNo : PSS-LIMS-SS-008

PHARMA SOFT

STABILITY MANAGEMENT

LIMS USER

Sample Acknowledgment

Action

Sample Acknowledgment Details

• Enter Quantity withdraw and select UOM from drop down.click on Select to select the Approver.

• If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish

• Enter E- Signature Password and Click On Finish.

Quantity Withdraw :

2

UOM :

Gm

Approval From:

Select

Used For:

Maximum allowed characters should be 100


Select an Action*:

☐ Acknowledge

☐ Terminate

I IMS 2.0

All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Sample Acknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

9. Enter Used For

SSNo : PSS-LIMS-SS-009

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Sample Acknowledgment

Action
Sample Acknowledgment Details

- Enter Quantity withdraw and select UOM from drop down. click on Select to select the Approver.
- If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish
- Enter E- Signature Password and Click On Finish.

Quantity Withdraw :
UOM :
Approval From: Select

Used For:

Select an Action*:
☐ Acknowledge
☐ Terminate

IMS 2.0
All © rights Reserved PharmaSoftSol

10. Click on Select

SSNo : PSS-LIMS-SS-010

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Sample Acknowledgment

Action
Sample Acknowledgment Details

- Enter Quantity withdraw and select UOM from drop down. click on Select to select the Approver.
- If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish
- Enter E- Signature Password and Click On Finish.

Quantity Withdraw :
UOM :
Approval From: Select

Used For:

Select an Action*:
☐ Acknowledge
☐ Terminate

Locations

organization

First Name:
Last Name:
Search Clear


Users

First Name:
Last Name:

Go to page: Row count:

Select Cancel

IMS 2.0
All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Sample Acknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

11. Click on Location

SSNo : PSS-LIMS-SS-011

PHARMA SOFT

STABILITY MANAGEMENT

LIMS USER

Sample Acknowledgment

Action

Sample Acknowledgment

• Enter Quantity

• If select Acknowledgment

• Enter E- Signature

Quantity

Withdrawal

Used For

Locations

organization

Plant-1

First Name:

Last Name:

Search

Clear

Users

First Name

Last Name

No data available!

Go to page: 1 Row count: 5

Select

Cancel

LIMS 2.0

All © rights Reserved PharmaSoftSol

12. Select a Record

SSNo : PSS-LIMS-SS-012

PHARMA SOFT
STABILITY MANAGEMENT

 LIMS USER

Sample Acknowledgement

Action Sample Acknowledgement

- Enter Quantity
- If select Acknowledgement
- Enter E-Signature

Quantity
Withdrawal:

Used For:

Locations


- organization
 - Plant-1
 - QA
 - QC
 - PF

First Name:
Last Name:

Users	
First Name	Last Name
LIMS	USER
LIMS	REVIEWER
LIMS	APPROVER

<< < 1 > >>
Go to page: Row count:
Showing 1-3 of 3

LIMS 2.0
All © rights Reserved PharmaSoftSol

	Test Case Name	Create Sample Acknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

13. Click on Select

SSNo : PSS-LIMS-SS-013

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Sample Acknowledgment

Action
Sample Acknowledgment Details

- Enter Quantity withdraw and select UOM from drop down. click on Select to select the Approver.
- If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish
- Enter E- Signature Password and Click On Finish.

Quantity Withdraw :
UOM :
Approval From: Select

Used For:

Select an Action*:
☐ Acknowledge
☐ Terminate

IMS 2.0
All © rights Reserved PharmaSoftSol

14. Click on Acknowledge

SSNo : PSS-LIMS-SS-014

PHARMA SOFT
STABILITY MANAGEMENT
LIMS USER

Sample Acknowledgment

Action
Sample Acknowledgment Details


- Enter Quantity withdraw and select UOM from drop down. click on Select to select the Approver.
- If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish
- Enter E- Signature Password and Click On Finish.

Quantity Withdraw :
UOM :
Approval From: Select

Used For:

Select an Action*:
☒ Acknowledge
☐ Terminate

IMS 2.0
All © rights Reserved PharmaSoftSol

	Test Case Name	Create Sample Acknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

15.Click on Finish

SSNo : PSS-LIMS-SS-015

PHARMA SOFT

STABILITY MANAGEMENT

LIMS USER

Sample Acknowledgment

Action

Sample Acknowledgment Details

• Enter Quantity withdraw and select UOM from drop down.click on Select to select the Approver.

• If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish

• Enter E- Signature Password and Click On Finish.

Quantity Withdraw :

2

Used For:

Sampling

Select an Action*:

Acknow

Validate E-Signature

E-Signature Password* :

Submit

Cancel

LIMS APPROVER

Select

LIMS 2.0

All © rights Reserved PharmaSoftSol

16.Enter E-Signature Password

SSNo : PSS-LIMS-SS-016

PHARMA SOFT

STABILITY MANAGEMENT

LIMS USER

Sample Acknowledgment

Action

Sample Acknowledgment Details

• Enter Quantity withdraw and select UOM from drop down.click on Select to select the Approver.

• If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish

• Enter E- Signature Password and Click On Finish.

Quantity Withdraw :

2

Used For:

Sampling

Select an Action*:

Acknow

Validate E-Signature

E-Signature Password* :

Submit


Cancel

LIMS APPROVER

Select

LIMS 2.0

All © rights Reserved PharmaSoftSol

	Test Case Name	Create Sample Acknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

17.Click on Submit

SSNo : PSS-LIMS-SS-017

PHARMA SOFT

STABILITY MANAGEMENT

LIMS USER

Sample Acknowledgment

Action

Sample Acknowledgment Details

• Enter Quantity withdraw and select UOM from drop down.click on Select to select the Approver.

• If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish

• Enter E- Signature Password and Click On Finish.

Quantity Withdraw :

2

Used For:

Sampling

Select an Action*:

Acknowledgment

Validate E-Signature

E-Signature Password* :

Saving...

Submit

Cancel

LIMS APPROVER

Select

LIMS 2.0

All © rights Reserved PharmaSoftSol

18.Click on OK button

SSNo : PSS-LIMS-SS-018

PHARMA SOFT

STABILITY MANAGEMENT

LIMS USER

Sample Acknowledgment

Action

Sample Acknowledgment Details

• Enter Quantity withdraw and select UOM from drop down.click on Select to select the Approver.

• If select Acknowledgment then enter Used for comments , or if select Terminate, then enter Comment if Terminate, click finish

• Enter E- Signature Password and Click On Finish.

Quantity Withdraw :

2

UOM :

Used For:

Sampling

Select an Action*:

Acknowledgment

Success

Submitted successfully.

OK

LIMS APPROVER


Select

LIMS 2.0

All © rights Reserved PharmaSoftSol

© SmartPharmaSoftSol

Page 9 of 10


	Test Case Name	Create Sample Acknowledgement
	Test Case Number	LIMS-SM-013
	Actual Result	Pass

19. Click on User Name

SSNo : PSS-LIMS-SS-019

PHARMA **SOFT SOL**
COMMITTED TO PAPERLESS COMPLIANCE

STABILITY MANAGEMENT

 **LIMS USER**

Sample Acknowledgment


Action
Sample Acknowledgment Details

. Select Sample Acknowledgment, Click On Search, Select The Record and Click on Next.
 . If Reverted Sample Acknowledgments by Approver is Required, Then Select Reverted Sample Acknowledgments by Approver.

☐ Sample Acknowledgment
 ☒ Reverted Sample Acknowledgments by Approver

S.No	Protocol Id	Product Name	Specification Number	AR Number	Protocol Number	Storage Conditions	Stability Conditions	Packing Type	Batch Number	Start Date	RevertedBy/On	Comments
1	7	Aspirin	SNSS 001	ASS220013	PS0003	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	Box	Product Batch No 1	11/06/2022	LIMS APPROVER 13/07/2022	View

<< < 1 > >>
 Go to page: Row count:
Showing 1-1 of 1



Logout

1 LIMS 2 0

All © rights Reserved PharmaSoftSol

20. Click on Logout

SSNo : PSS-LIMS-SS-020

PHARMA **SOFT SOL**
COMMITTED TO PAPERLESS COMPLIANCE

LIMS

User Name

Password

Module

Select

Login

Reset