	Test Case Name	Approve Protocol
	Test Case Number	LIMS-SM-005
	Actual Result	Pass

1. Enter the username, password, Select Module and click on login button


SSNo : PSS-LIMS-SS-001



The login form is titled "LIMS" and features the PharmaSoftSol logo. It contains three input fields: "User Name" with the value "LIMS_APPROVER", "Password" with masked characters "*****", and "Module" with a dropdown menu set to "STABILITY MANAGEMENT". Below the fields are "Login" and "Reset" buttons.

2. Click on Protocol

SSNo : PSS-LIMS-SS-002

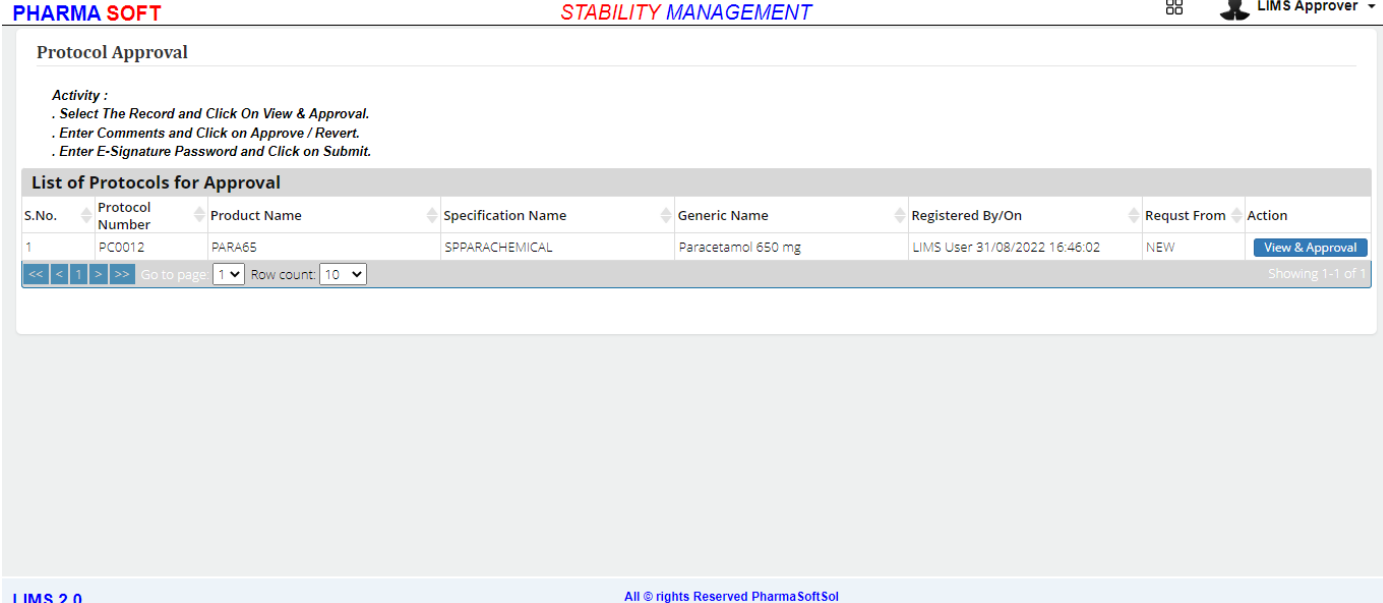


The page header shows "PHARMA SOFT" and "STABILITY MANAGEMENT". The user is logged in as "LIMS Approver". The main section is titled "Protocol Approval" and includes instructions: "Activity : . Select The Record and Click On View & Approval. . Enter Comments and Click on Approve / Revert. . Enter E-Signature Password and Click on Submit." Below this is a table titled "List of Protocols for Approval".

S.No.	Protocol Number	Product Name	Specification Name	Generic Name	Registered By/On	Request From	Action
1	PC0012	PARA65	SPPARACHEMICAL	Paracetamol 650 mg	LIMS User 31/08/2022 16:46:02	NEW	View & Approval

Navigation: << < 1 > >> Go to page: 1 Row count: 10 Showing 1-1 of 1

Footer: LIMS 2.0 All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Approve Protocol
	Test Case Number	LIMS-SM-005
	Actual Result	Pass

3. Select a Record And Click on View & Approval

SSNo : PSS-LIMS-SS-003

PHARMA SOFT
STABILITY MANAGEMENT
88
LIMS Approver

rotocol Stability Plan

Procedure :

↶

↷

Formats

B

I

≡

≡

≡

≡

≡

≡

≡

≡

rotocol Procedure

Location:

Plant-1

Specification Number:

SPPARACH001

Generic Name:

Paracetamol 650 mg

Test Planner Version:

0

Product Name:

PARA65

Specification Name:

SPPARACHEMICAL

STP Number:

STPCHPARA001

Customer Name:

BA

4. Enter Comments

SSNo : PSS-LIMS-SS-004

PHARMA SOFT
STABILITY MANAGEMENT
88
LIMS Approver


Packing type	No	Batch size	Manufactured Date	Expiry Date	Storage Condition	Orientation	Stability Condition	Type
Blister	44	100	08/31/2022	08/31/2027	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	MONTH

Comments*

Approved

Approve

Revert


 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Approve Protocol
	Test Case Number	LIMS-SM-005
	Actual Result	Pass

5. Click on Approve

SSNo : PSS-LIMS-SS-005

PHARMA SOFT

STABILITY MANAGEMENT

 LIMS Approver

Packing type	No	Batch size	Manufactured Date	Expiry Date	Storage Condition	Orientation	Stability Condition	Type
Blister	44	100	08/31/2022	08/31/2027	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	MONTH

Comments*:

Approved

Validate E-Signature

E-Signature Password*:

Submit

Cancel

Approve


Revert

6. Enter E-Signature Password

SSNo : PSS-LIMS-SS-006

PHARMA SOFT

STABILITY MANAGEMENT

 LIMS Approver

Packing type	No	Batch size	Manufactured Date	Expiry Date	Storage Condition	Orientation	Stability Condition	Type
Blister	44	100	08/31/2022	08/31/2027	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	MONTH

Comments*:

Approved

Validate E-Signature


E-Signature Password*:

Submit

Cancel

Approve

Revert

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Approve Protocol
	Test Case Number	LIMS-SM-005
	Actual Result	Pass

7. Click on Submit

SSNo : PSS-LIMS-SS-007

PHARMA SOFT
STABILITY MANAGEMENT
☰
👤 LIMS Approver

Packing type	No	Batch size	Manufactured Date	Expiry Date	Storage Condition	Orientation	Stability Condition	Type
Blister	44	100	08/31/2022	08/31/2027	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	MONTH

Comments*:

Approved

Validate E-Signature

E-Signature Password*:

Saving...

Submit

Cancel

Approve

Revert

8. Click on OK button

SSNo : PSS-LIMS-SS-008

PHARMA SOFT
STABILITY MANAGEMENT
☰
👤 LIMS Approver

Packing Type	Product Batch No	Batch Size	Manufactured Date	Expiry Date	Storage Condition	Orientation	Stability Condition	Interval Type	Int
Blister	44	100	08/31/2022	08/31/2027	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	MONTH	

Comments*:

Approved


Success

Submitted successfully.

OK

Approve


Revert

	Test Case Name	Approve Protocol
	Test Case Number	LIMS-SM-005
	Actual Result	Pass

9. Click on User Name


SSNo : PSS-LIMS-SS-009

PHARMA SOFT SOL
STABILITY MANAGEMENT

 LIMS Approver

Protocol Approval

Activity :
 . Select The Record and Click On View & Approval.
 . Enter Comments and Click on Approve / Revert.
 . Enter E-Signature Password and Click on Submit.



List of Protocols for Approval

S.No.	Protocol Number	Product Name	Specification Name	Generic Name	Registered By/On	Request From	Action
No data available!							

Row count: 10

LIMS 2.0
 All © rights Reserved PharmaSoftSol

10. Click on Logout

SSNo : PSS-LIMS-SS-010

LIMS



User Name

Password

Module

Select

Login

Reset