	Test Case Name	Reject Protocol
	Test Case Number	PSS-LIMS-005
	Actual Result	Pass

1. Enter the username, password, Select Module and click on login button


SSNo : PSS-LIMS-SS-001



The screenshot shows the LIMS login interface. At the top, the 'LIMS' title is displayed in yellow. Below it is the PharmaSoftSol logo. The login form contains three input fields: 'User Name' with the value 'LIMS_APPROVER', 'Password' with masked characters '*****', and 'Module' with a dropdown menu set to 'STABILITY MANAGEMENT'. At the bottom of the form are two buttons: 'Login' and 'Reset'.

2. Click on Protocol

SSNo : PSS-LIMS-SS-002

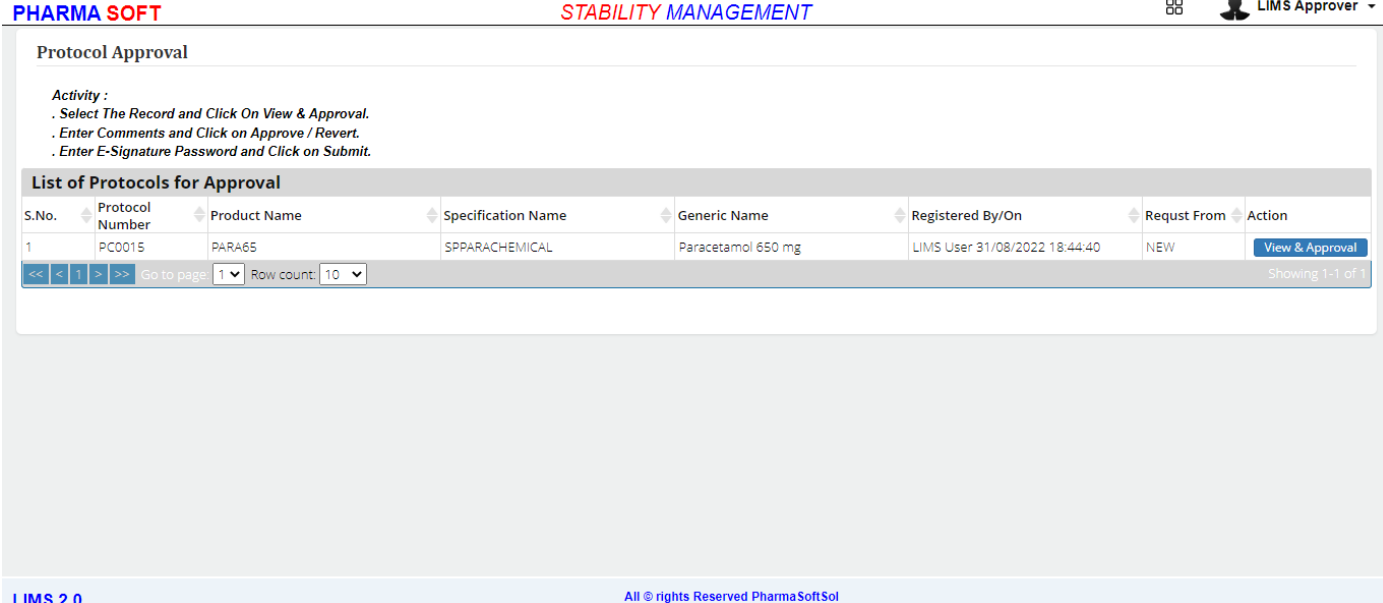


The screenshot displays the 'Protocol Approval' page within the LIMS system. The header includes the PharmaSoftSol logo, the title 'STABILITY MANAGEMENT', and a user profile icon labeled 'LIMS Approver'. The main content area is titled 'Protocol Approval' and includes instructions for activity: 'Select The Record and Click On View & Approval.', 'Enter Comments and Click on Approve / Revert.', and 'Enter E-Signature Password and Click on Submit.'. Below the instructions is a table titled 'List of Protocols for Approval'.

S.No.	Protocol Number	Product Name	Specification Name	Generic Name	Registered By/On	Request From	Action
1	PC001S	PARA6S	SPPARACHEMICAL	Paracetamol 650 mg	LIMS User 31/08/2022 18:44:40	NEW	View & Approval

Below the table, there are navigation controls: '<<' '<' '1' '>' '>>' 'Go to page: 1' 'Row count: 10' and 'Showing 1-1 of 1'.

The footer of the page contains the text 'LIMS 2.0' and 'All © rights Reserved PharmaSoftSol'.

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Reject Protocol
	Test Case Number	PSS-LIMS-005
	Actual Result	Pass

3. Select a Record And Click on View & Approval

SSNo : PSS-LIMS-SS-003

PHARMA SOFT
STABILITY MANAGEMENT
88
LIMS Approver

rotocol Stability Plan

Procedure :

↶

↷

Formats

B

I

≡

≡

≡

≡

≡

≡

≡

≡

rotocol Procedure

Location:

Plant-1

Specification Number:

SPPARACH001

Generic Name:

Paracetamol 650 mg

Test Planner Version:

0

Product Name:

PARA65

Specification Name:

SPPARACHEMICAL

STP Number:

STPCHPARA001

Customer Name:

BA

4. Enter Comments

SSNo : PSS-LIMS-SS-004

PHARMA SOFT
STABILITY MANAGEMENT
88
LIMS Approver


Packing type	No	Batch size	Manufactured Date	Expiry Date	Storage Condition	Orientation	Stability Condition	Type
Blister	44	100	08/31/2022	08/31/2027	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	MONTH

Comments*:

Rejected

Approve

Revert


 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Reject Protocol
	Test Case Number	PSS-LIMS-005
	Actual Result	Pass

5. Click on Reject

SSNo : PSS-LIMS-SS-005

PHARMA SOFT

STABILITY MANAGEMENT

 LIMS Approver

Packing type	No	Batch size	Manufactured Date	Expiry Date	Storage Condition	Orientation	Stability Condition	Type
Blister	44	100	08/31/2022	08/31/2027	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	MONTH

Comments*:

Rejected

Validate E-Signature

E-Signature Password*:

Submit

Cancel

Approve


Revert

6. Enter E-Signature Password

SSNo : PSS-LIMS-SS-006

PHARMA SOFT

STABILITY MANAGEMENT

 LIMS Approver

Packing type	No	Batch size	Manufactured Date	Expiry Date	Storage Condition	Orientation	Stability Condition	Type
Blister	44	100	08/31/2022	08/31/2027	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	MONTH

Comments*:

Rejected

Validate E-Signature


E-Signature Password*:

Submit

Cancel

Approve

Revert

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Reject Protocol
	Test Case Number	PSS-LIMS-005
	Actual Result	Pass

7.Click on Submit

SSNo : PSS-LIMS-SS-007

PHARMA SOFT
STABILITY MANAGEMENT
☰
👤 LIMS Approver ▾

Packing type	No	Batch size	Manufactured Date	Expiry Date	Storage Condition	Orientation	Stability Condition	Type
Blister	44	100	08/31/2022	08/31/2027	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	MONTH

Comments*:

Rejected

Approve
Revert

8.Click on OK button

SSNo : PSS-LIMS-SS-008

PHARMA SOFT
STABILITY MANAGEMENT
☰
👤 LIMS Approver ▾


Packing Type	Product Batch No	Batch Size	Manufactured Date	Expiry Date	Storage Condition	Orientation	Stability Condition	Interval Type	Int
Blister	44	100	08/31/2022	08/31/2027	40°C ± 2°C / 75%RH ± 5%RH	UP RIGHT	ACCELERATED	MONTH	

Comments*:

Rejected

Success
Submitted successfully.
OK


Approve
Revert

	Test Case Name	Reject Protocol
	Test Case Number	PSS-LIMS-005
	Actual Result	Pass

9.Click on User Name


SSNo : PSS-LIMS-SS-009

PHARMA SOFT SOL
STABILITY MANAGEMENT

 LIMS Approver

Protocol Approval

Activity :
 . Select The Record and Click On View & Approval.
 . Enter Comments and Click on Approve / Revert.
 . Enter E-Signature Password and Click on Submit.



List of Protocols for Approval

S.No.	Protocol Number	Product Name	Specification Name	Generic Name	Registered By/On	Request From	Action
No data available!							

Row count: 10

LIMS 2.0
 All © rights Reserved PharmaSoftSol

10.Click on Logout

SSNo : PSS-LIMS-SS-010

LIMS



User Name

Password

Module