
	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

1. Enter the username, password, Select Module and click on login button


SSNo : PSS-LIMS-SS-001

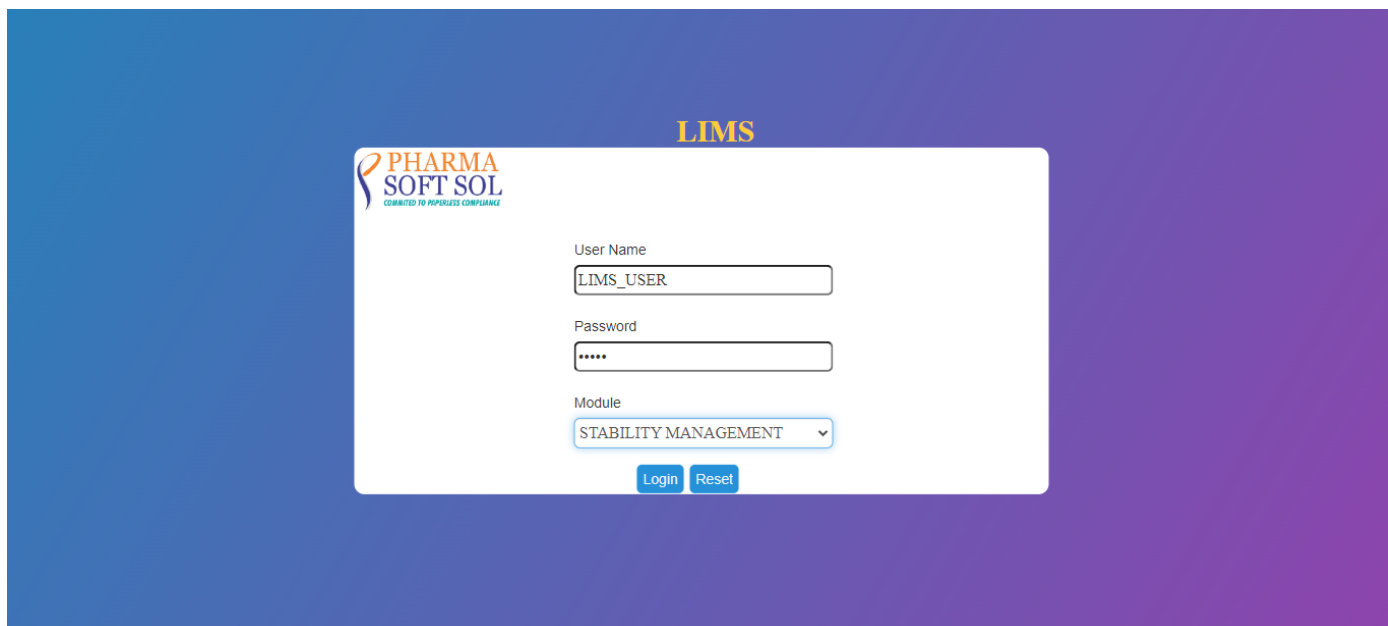

The screenshot shows the LIMS login interface. At the top, the 'LIMS' title is displayed in yellow. Below it, the PharmaSoftSol logo is visible. The login form contains three input fields: 'User Name' with the value 'LIMS_USER', 'Password' with masked characters '*****', and 'Module' with a dropdown menu set to 'STABILITY MANAGEMENT'. At the bottom of the form are two buttons: 'Login' and 'Reset'.

1. Click on Loading

SSNo : PSS-LIMS-SS-001

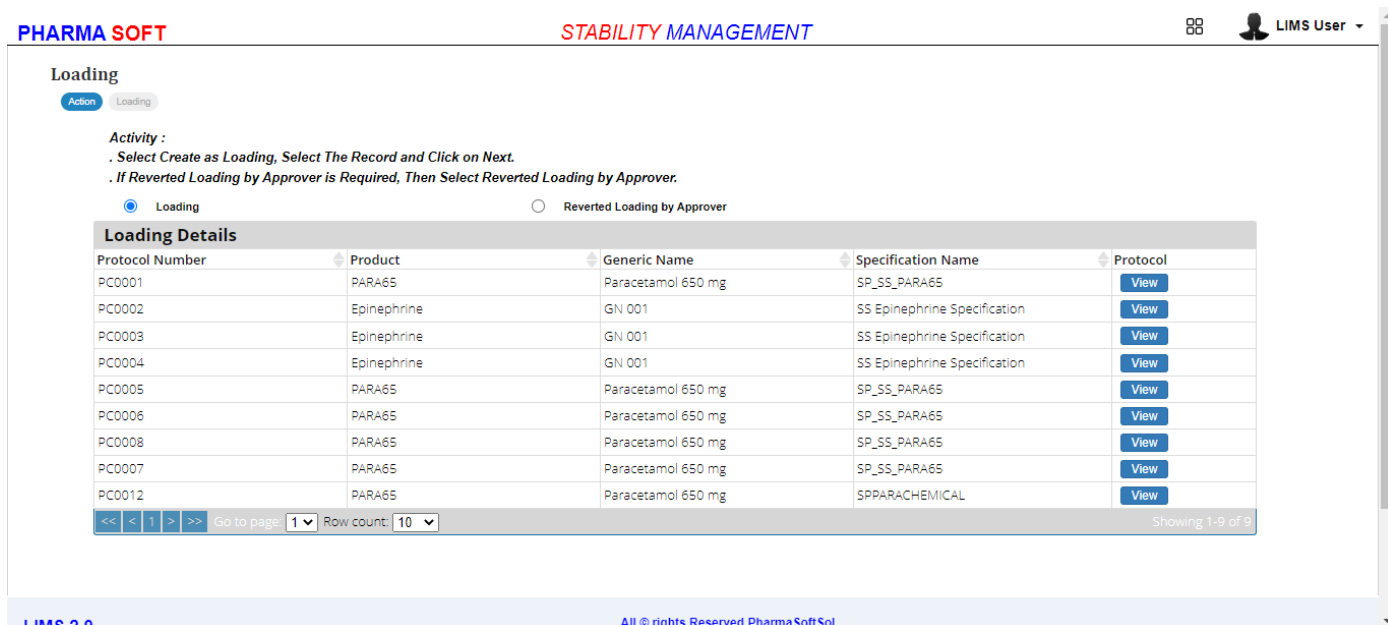

The screenshot displays the 'Loading' section of the LIMS application. The header includes the PharmaSoftSol logo, the title 'STABILITY MANAGEMENT', and a user profile icon labeled 'LIMS User'. The 'Loading' section has two tabs: 'Action' (selected) and 'Loading'. Below the tabs, an 'Activity' section provides instructions: 'Select Create as Loading, Select The Record and Click on Next. If Reverted Loading by Approver is Required, Then Select Reverted Loading by Approver.' There are two radio buttons: 'Loading' (selected) and 'Reverted Loading by Approver'. The main area is a table titled 'Loading Details' with columns: Protocol Number, Product, Generic Name, Specification Name, and Protocol. The table contains 12 rows of data, each with a 'View' button in the Protocol column. At the bottom, there is a pagination bar with navigation arrows, 'Go to page: 1', 'Row count: 10', and 'Showing 1-9 of 9'.

Protocol Number	Product	Generic Name	Specification Name	Protocol
PC0001	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0002	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0003	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0004	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0005	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0006	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0008	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0007	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0012	PARA65	Paracetamol 650 mg	SPPARACHEMICAL	View

2. Select a Record

SSNo : PSS-LIMS-SS-002

Test Case Name	Create Loading
Test Case Number	LIMS-SM-011
Actual Result	Pass

PHARMA SOFT

STABILITY MANAGEMENT

LIMS User

Loading

Action

Loading

Activity :

. Select Create as Loading, Select The Record and Click on Next.

. If Reverted Loading by Approver is Required, Then Select Reverted Loading by Approver.

Loading

Reverted Loading by Approver

Loading Details

Protocol Number	Product	Generic Name	Specification Name	Protocol
PC0001	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0002	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0003	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0004	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0005	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0006	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0008	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0007	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0012	PARA65	Paracetamol 650 mg	SPPARACHEMICAL	View

<<

<

1

>

>>

Go to page

1

Row count:

10


Showing 1-9 of 9

IMS 2.0

All © rights Reserved PharmaSoftSol

© SmartPharmaSoftSol

Page 2 of 13

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

3. Click on Next

SSNo : PSS-LIMS-SS-003

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading

Action Loading

ocol

EW

EW

EW

EW

EW

EW

EW

EW

EW

Showing 1-9 of 9

IMS 2.0
All © rights Reserved PharmaSoftSol

4. Enter Analysis Quantity

SSNo : PSS-LIMS-SS-004

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading

Action Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.

. Select Approval User From Click On Select and Click On Finish.

. Enter E-Signature Password and Click On Submit.


Loading Details

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2		Selec			Sele	Evaluate		Select

Approval From*:

Select

IMS 2.0
All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

5. Enter Reference Quantity

SSNo : PSS-LIMS-SS-005

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading

Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
. Select Approval User From Click On Select and Click On Finish.
. Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Sele			Sele	Evaluate		Select

Approval From*:
Select

IMS 2.0
All © rights Reserved PharmaSoftSol

6. Select UOM

SSNo : PSS-LIMS-SS-006

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading


Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
. Select Approval User From Click On Select and Click On Finish.
. Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm			Sele	Evaluate		Select

Approval From*:
Select

IMS 2.0
All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

7. Enter Analysis Quantity Count

SSNo : PSS-LIMS-SS-007

PHARMA SOFT
STABILITY MANAGEMENT
☰
LIMS User ▾

Loading

Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
. Select Approval User From Click On Select and Click On Finish.
. Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm ▾	3		Sel ▾	Evaluate		Select

Approval From*: Select

IMS 2.0
All © rights Reserved PharmaSoftSol

8. Enter Reference Quantity Count

SSNo : PSS-LIMS-SS-008

PHARMA SOFT
STABILITY MANAGEMENT
☰
LIMS User ▾

Loading


Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
. Select Approval User From Click On Select and Click On Finish.
. Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm ▾	3	2	Sel ▾	Evaluate		Select

Approval From*: Select

IMS 2.0
All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

9. Select UOM

SSNo : PSS-LIMS-SS-009

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading

Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
 . Select Approval User From Click On Select and Click On Finish.
 . Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	<input type="text" value="2"/>	<input type="text" value="3"/>	Gm	<input type="text" value="3"/>	<input type="text" value="2"/>	Gm	<input type="button" value="Evaluate"/>	<input type="text"/>	<input type="button" value="Select"/>

Approval From*:

IIMS 2.0
 All © rights Reserved PharmaSoftSol

10. Click on Start Date

SSNo : PSS-LIMS-SS-010

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading


Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
 . Select Approval User From Click On Select and Click On Finish.
 . Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	<input type="text" value="2"/>	<input type="text" value="3"/>	Gm	<input type="text" value="3"/>	<input type="text" value="2"/>	Gm	<input type="button" value="Evaluate"/>	<input type="text" value="31-08-2022"/>	<input type="button" value="Select"/>


Approval From*:

IIMS 2.0
 All © rights Reserved PharmaSoftSol

	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

11.Click on Select

SSNo : PSS-LIMS-SS-011

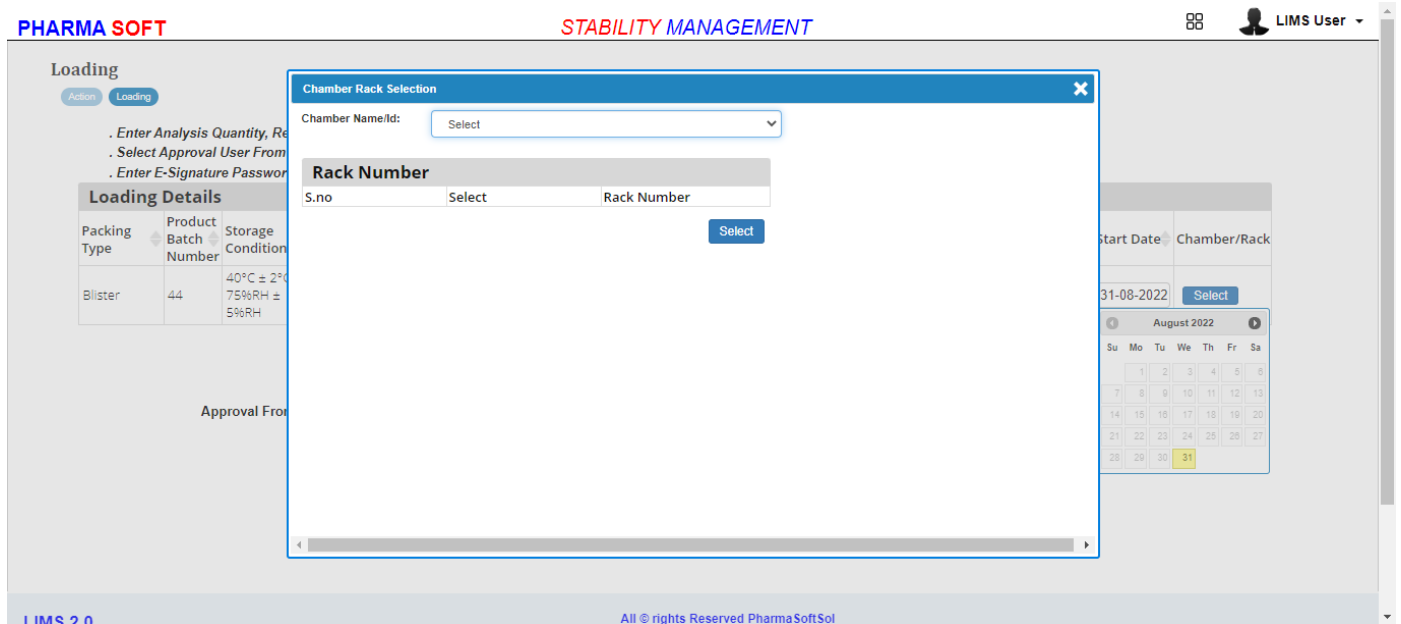

PHARMA SOFT **STABILITY MANAGEMENT** LIMS User

Loading

Enter Analysis Quantity, Re
Select Approval User From
Enter E-Signature Passwor

Loading Details

Packing Type	Product Batch Number	Storage Condition
Blister	44	40°C ± 2°C 75%RH ± 5%RH

Approval Fro

Chamber Rack Selection

Chamber Name/Id: Select

Rack Number

S.no	Select	Rack Number
------	--------	-------------

Select

Start Date: 31-08-2022 Chamber/Rack: Select


August 2022

Su	Mo	Tu	We	Th	Fr	Sa
1	2	3	4	5	6	
7	8	9	10	11	12	13
14	15	16	17	18	19	20
21	22	23	24	25	26	27
28	29	30	31			

LIMS 2.0 All © rights Reserved PharmaSoftSol

12.Select Chamber

SSNo : PSS-LIMS-SS-012

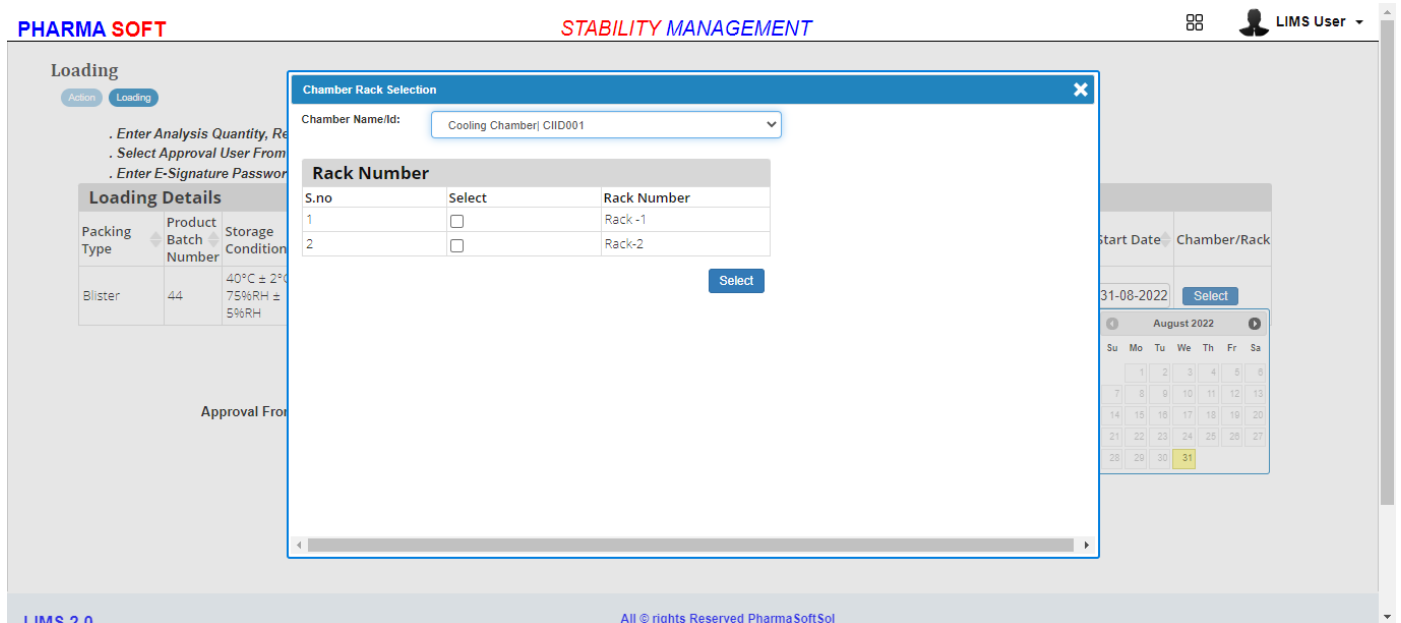

PHARMA SOFT **STABILITY MANAGEMENT** LIMS User

Loading

Enter Analysis Quantity, Re
Select Approval User From
Enter E-Signature Passwor

Loading Details

Packing Type	Product Batch Number	Storage Condition
Blister	44	40°C ± 2°C 75%RH ± 5%RH

Approval Fro

Chamber Rack Selection

Chamber Name/Id: Cooling Chamber| CIID001

Rack Number

S.no	Select	Rack Number
1	<input type="checkbox"/>	Rack -1
2	<input type="checkbox"/>	Rack-2


Select

Start Date: 31-08-2022 Chamber/Rack: Select

August 2022

Su	Mo	Tu	We	Th	Fr	Sa
1	2	3	4	5	6	
7	8	9	10	11	12	13
14	15	16	17	18	19	20
21	22	23	24	25	26	27
28	29	30	31			

LIMS 2.0 All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

13. Select Check Box

SSNo : PSS-LIMS-SS-013

PHARMA SOFT **STABILITY MANAGEMENT** LIMS User

Loading
 Action Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
 . Select Approval User From Click On Select and Click On Finish.
 . Enter E-Signature Password and Click On Submit.

Loading Details

Packing Type	Product Batch Number	Storage Conditions
Blister	44	40°C ± 2°C / 75%RH ± 5%RH

Approval From:

Chamber Rack Selection

Chamber Name/Id:

S.no	Select	Rack Number
1	<input checked="" type="checkbox"/>	Rack-1
2	<input type="checkbox"/>	Rack-2

Start Date:

August 2022

Su	Mo	Tu	We	Th	Fr	Sa
1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30	31				

IMS 2.0 All © rights Reserved PharmaSoftSol

14. Click on Select

SSNo : PSS-LIMS-SS-014

PHARMA SOFT **STABILITY MANAGEMENT** LIMS User

Loading
 Action Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
 . Select Approval User From Click On Select and Click On Finish.
 . Enter E-Signature Password and Click On Submit.

Loading Details


Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	<input type="text" value="2"/>	<input type="text" value="3"/>	Gm	<input type="text" value="3"/>	<input type="text" value="2"/>	Gm	<input type="button" value="Evaluate"/>	<input type="text" value="31-08-2022"/>	<input type="button" value="Select"/>

Approval From*:

August 2022

Su	Mo	Tu	We	Th	Fr	Sa
1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30	31				

IMS 2.0 All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

15. Click on Select

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading

Action Loading

- . Enter Anal
- . Select App
- . Enter E-Sig

Loading De

Packing Type	Pro Bat Nu
Blister	44

Locations

organization

First Name:

Last Name:

Search

Clear

Users

First Name	Last Name
No data available!	

Go to page: 1
Row count: 5

Select

Cancel

16. Click on Location

PHARMA SOFT

STABILITY MANAGEMENT

Loading

Action Loading

. Enter Anal

. Select App

. Enter E-Sig

Loading De

Packing Type

Blister

44

Locations

organization

Plant-1

QA

QC

EN

First Name:

Last Name:

Search Clear

Users

First Name

Last Name

No data available!

Go to page: 1 Row count: 5

Select Cancel

select

22

Th Fr Sa

4 5 6


11 12 13

18 19 20

25 26 27

LIMS 2.0

All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

17. Select a Record

SSNo : PSS-LIMS-SS-017

PHARMA SOFT **STABILITY MANAGEMENT** LIMS User

Loading

Action Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.

. Select Approval User From Click On Select and Click On Finish.

. Enter E-Signature Password and Click On Submit.

Loading Details

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm	3	2	Gm		31-08-2022	

Approval From*: Select

Locations

organization

Plant-1

QA

QC

EN

First Name: Last Name:

Search Clear

Users

First Name	Last Name
LIMS	User
LIMS	Approver
LIMS	Reviewer
User	01
User	02

<< < 1 2 3 > >> Go to page: Row count: Showing 1-5 of 12

Select Cancel

IIMS 2.0 All © rights Reserved PharmaSoftSol

18. Click on Select

SSNo : PSS-LIMS-SS-018

PHARMA SOFT **STABILITY MANAGEMENT** LIMS User

Loading

Action Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.

. Select Approval User From Click On Select and Click On Finish.


. Enter E-Signature Password and Click On Submit.

Loading Details

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm	3	2	Gm		31-08-2022	

Approval From*: Select

IIMS 2.0 All © rights Reserved PharmaSoftSol

	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

21.Click on Submit

SSNo : PSS-LIMS-SS-021

PHARMA SOFT **STABILITY MANAGEMENT** ☰ LIMS User

Enter E-Signature Password and Click on Submit

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm	3	2	Gm	Evaluate	31-08-2022	Select

Approval From*: LIMS Approver

Validate E-Signature

E-Signature Password*:

Saving...

Submit Cancel

Previous Finish

LIMS 2.0 All © rights Reserved PharmaSoftSol

22.Click on OK button

SSNo : PSS-LIMS-SS-022

PHARMA SOFT **STABILITY MANAGEMENT** ☰ LIMS User

Enter E-Signature Password and Click on Submit

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm	3	2	Gm	Evaluate	31-08-2022	Select

Approval From*: LIMS Approver


Success

Submitted successfully.

OK


Previous Finish

LIMS 2.0 All © rights Reserved PharmaSoftSol


	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

23.Click on User Name

SSNo : PSS-LIMS-SS-023



STABILITY MANAGEMENT


 LIMS User

☒ Loading


☐ Reverted Loading by Approver

Protocol Number	Product	Generic Name	Specification Name	Protocol
PC0001	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0002	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0003	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0004	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0005	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0006	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0008	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0007	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View

[Previous](#)
[Next](#)



Logout



Go to page: Row count:


Showing 1-8 of 8

LIMS 2.0

All © rights Reserved PharmaSoftSol

24.Click on Logout

SSNo : PSS-LIMS-SS-024



LIMS

User Name

Password

Module

Login
Reset