
	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

1. Enter the username, password, Select Module and click on login button


SSNo : PSS-LIMS-SS-001

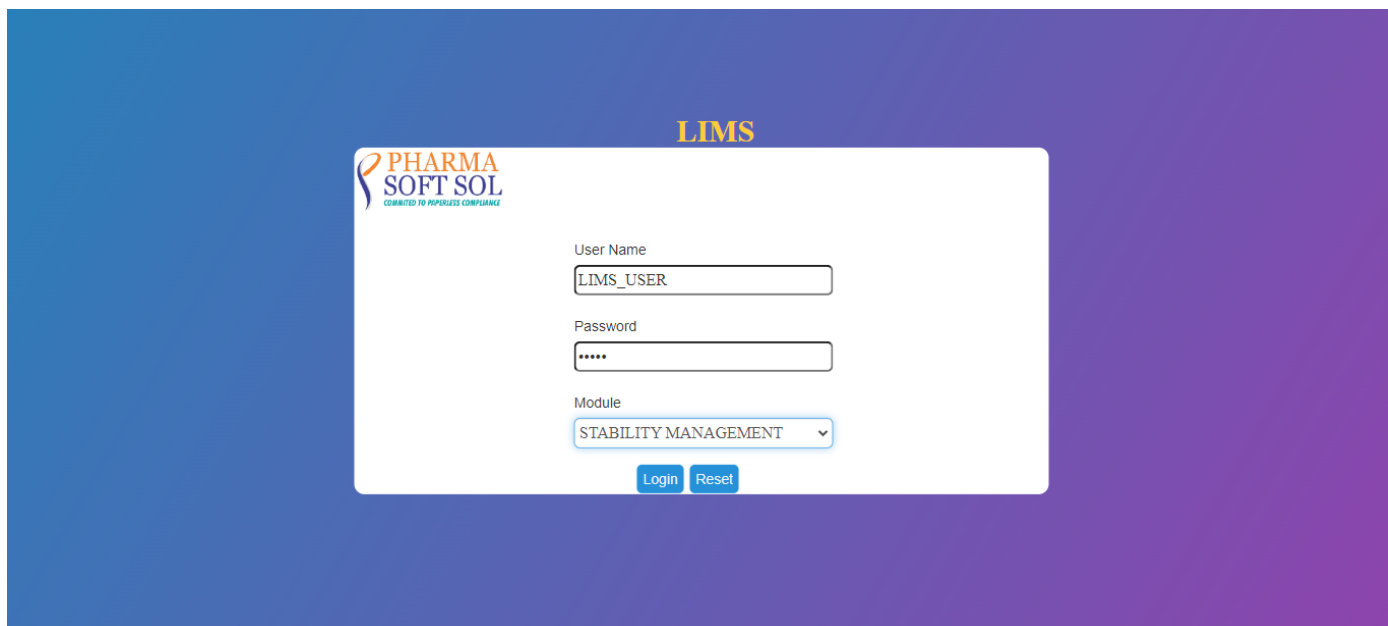

The screenshot shows the LIMS login interface. At the top, the 'LIMS' title is displayed in yellow. Below it, the PharmaSoftSol logo is visible. The login form contains three input fields: 'User Name' with the value 'LIMS_USER', 'Password' with masked characters '*****', and 'Module' with a dropdown menu set to 'STABILITY MANAGEMENT'. At the bottom of the form are two buttons: 'Login' and 'Reset'.

1. Click on Loading

SSNo : PSS-LIMS-SS-001

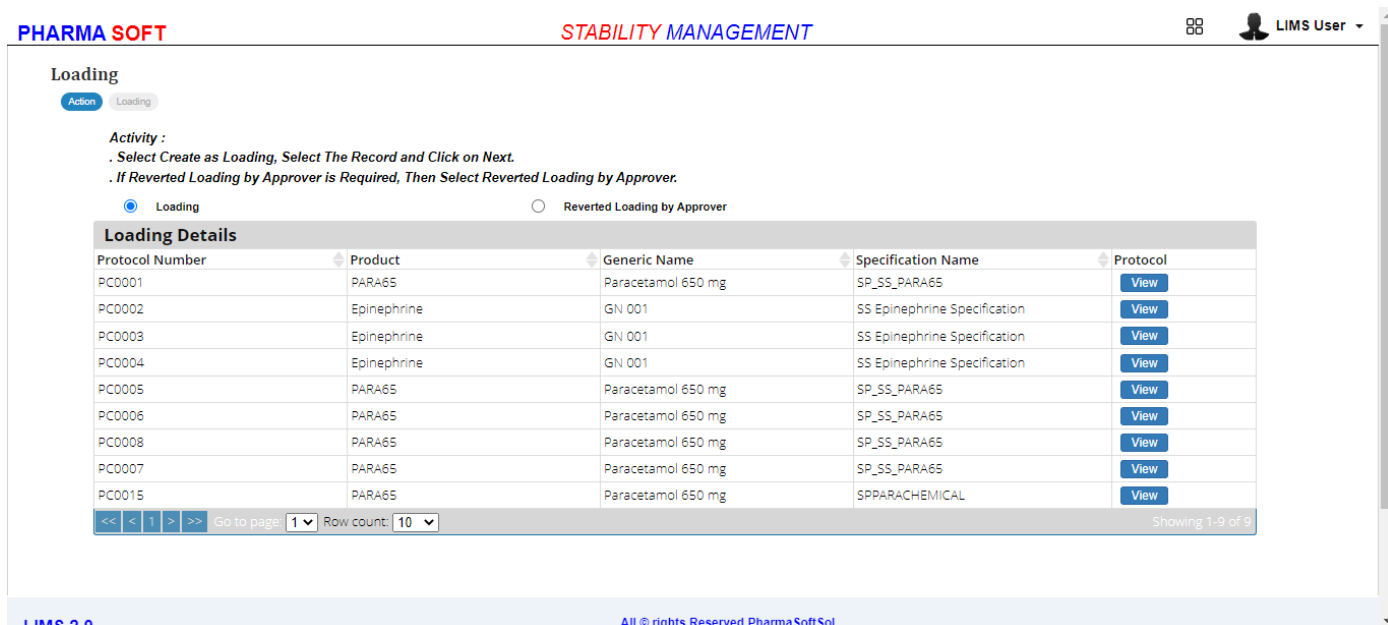

The screenshot displays the 'Loading' section of the LIMS application. The header includes the PharmaSoftSol logo, the title 'STABILITY MANAGEMENT', and a user profile icon labeled 'LIMS User'. The main content area is titled 'Loading' and includes an 'Action' button. Below this, there is an 'Activity' section with instructions: 'Select Create as Loading, Select The Record and Click on Next.' and 'If Reverted Loading by Approver is Required, Then Select Reverted Loading by Approver.' Two radio buttons are present: 'Loading' (selected) and 'Reverted Loading by Approver'. The 'Loading Details' table lists various records with columns for Protocol Number, Product, Generic Name, Specification Name, and Protocol. Each row has a 'View' button. The table data is as follows:

Protocol Number	Product	Generic Name	Specification Name	Protocol
PC0001	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0002	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0003	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0004	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0005	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0006	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0008	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0007	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0015	PARA65	Paracetamol 650 mg	SPPARACHEMICAL	View

At the bottom of the table, there is a pagination bar with navigation arrows, 'Go to page: 1', 'Row count: 10', and 'Showing 1-9 of 9'.

2. Select a Record

SSNo : PSS-LIMS-SS-002

Test Case Name	Create Loading
Test Case Number	LIMS-SM-011
Actual Result	Pass

PHARMA SOFT

STABILITY MANAGEMENT

LIMS User

Loading

Action

Loading

Activity :

. Select Create as Loading, Select The Record and Click on Next.

. If Reverted Loading by Approver is Required, Then Select Reverted Loading by Approver.

Loading

Reverted Loading by Approver

Loading Details

Protocol Number	Product	Generic Name	Specification Name	Protocol
PC0001	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0002	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0003	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0004	Epinephrine	GN 001	SS Epinephrine Specification	View
PC0005	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0006	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0008	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0007	PARA65	Paracetamol 650 mg	SP_SS_PARA65	View
PC0015	PARA65	Paracetamol 650 mg	SPPARACHEMICAL	View

<<

<

1

>

>>

Go to page

1

>

Row count:

10

>


Showing 1-9 of 9

IMS 2.0

All © rights Reserved PharmaSoftSol

© SmartPharmaSoftSol

Page 2 of 13

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

3. Click on Next

SSNo : PSS-LIMS-SS-003

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading

Action Loading

ation Name	Protocol
RA65	View
phrine Specification	View
phrine Specification	View
phrine Specification	View
RA65	View
RA65	View
RA65	View
RA65	View
RA65	View
CHEMICAL	View

Showing 1-9 of 9

IMS 2.0
All © rights Reserved PharmaSoftSol

4. Enter Analysis Quantity

SSNo : PSS-LIMS-SS-004

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading


Action Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
. Select Approval User From Click On Select and Click On Finish.
. Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2		Selec			Sele	Evaluate		Select

Approval From*:
[Select](#)

IMS 2.0
All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

5. Enter Reference Quantity

SSNo : PSS-LIMS-SS-005

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading

Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
. Select Approval User From Click On Select and Click On Finish.
. Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Sele			Sele	Evaluate		Select

Approval From*:
Select

IMS 2.0
All © rights Reserved PharmaSoftSol

6. Select UOM

SSNo : PSS-LIMS-SS-006

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading


Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
. Select Approval User From Click On Select and Click On Finish.
. Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm			Sele	Evaluate		Select

Approval From*:
Select

IMS 2.0
All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

7. Enter Analysis Quantity Count

SSNo : PSS-LIMS-SS-007

PHARMA SOFT
STABILITY MANAGEMENT
☰
LIMS User ▾

Loading
Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
 . Select Approval User From Click On Select and Click On Finish.
 . Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm ▾	3		Sel ▾	Evaluate		Select

Approval From*: Select

IMS 2.0
All © rights Reserved PharmaSoftSol

8. Enter Reference Quantity Count

SSNo : PSS-LIMS-SS-008

PHARMA SOFT
STABILITY MANAGEMENT
☰
LIMS User ▾


Loading
Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
 . Select Approval User From Click On Select and Click On Finish.
 . Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm ▾	3	2	Sel ▾	Evaluate		Select

Approval From*: Select

IMS 2.0
All © rights Reserved PharmaSoftSol

	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

9. Select UOM

SSNo : PSS-LIMS-SS-009

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading

Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
. Select Approval User From Click On Select and Click On Finish.
. Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm	3	2	Gm	Evaluate		Select

Approval From*: Select

IMS 2.0
All © rights Reserved PharmaSoftSol

10. Click on Start Date

SSNo : PSS-LIMS-SS-010

PHARMA SOFT
STABILITY MANAGEMENT
LIMS User

Loading


Action
Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
. Select Approval User From Click On Select and Click On Finish.
. Enter E-Signature Password and Click On Submit.

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm	3	2	Gm	Evaluate	31-08-2022	Select


Approval From*: Select

IMS 2.0
All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

11. Click on Select


SSNo : PSS-LIMS-SS-011

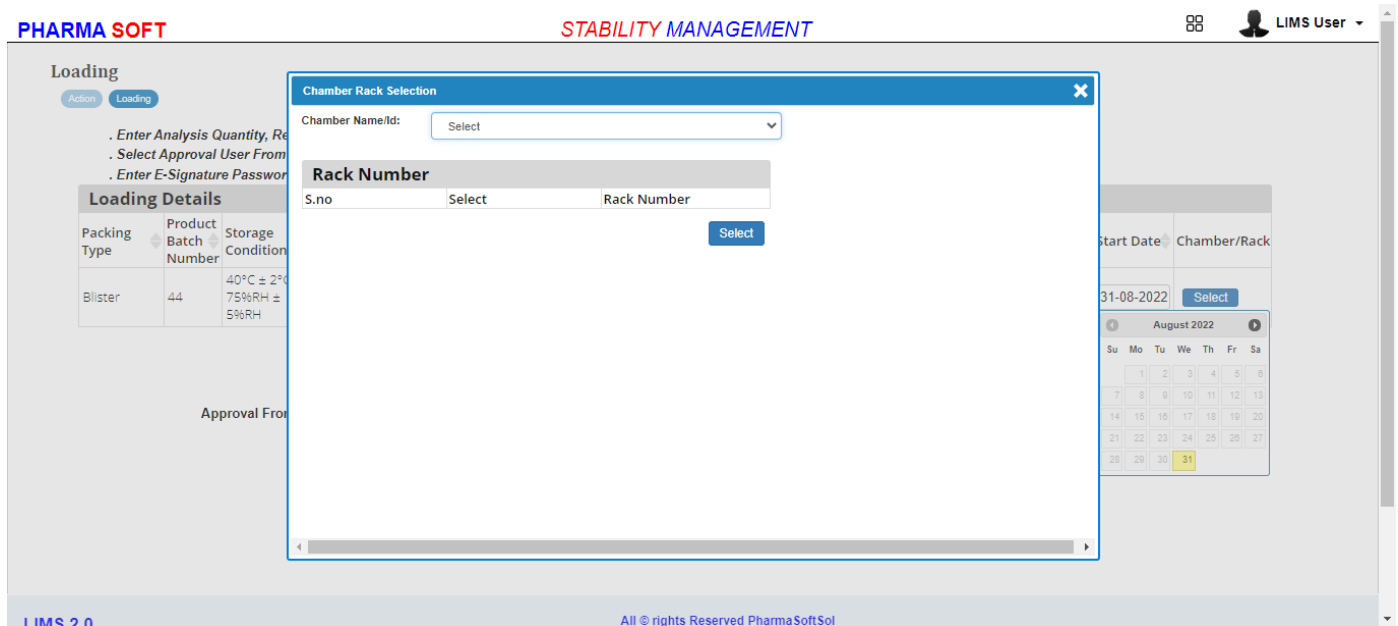


The screenshot shows the PHARMA SOFT STABILITY MANAGEMENT application. The 'Loading' section is active, displaying a table with 'Loading Details' including Packing Type (Blister), Product Batch Number (44), and Storage Condition (40°C ± 2°C, 75%RH ± 5%RH). A 'Chamber Rack Selection' dialog box is open, showing a 'Chamber Name/Id' dropdown set to 'Select' and a 'Rack Number' table with columns 'S.no' and 'Rack Number'. A 'Select' button is visible in the dialog. In the background, a calendar for August 2022 is visible, with the date 31-08-2022 selected.

12. Select Chamber

SSNo : PSS-LIMS-SS-012

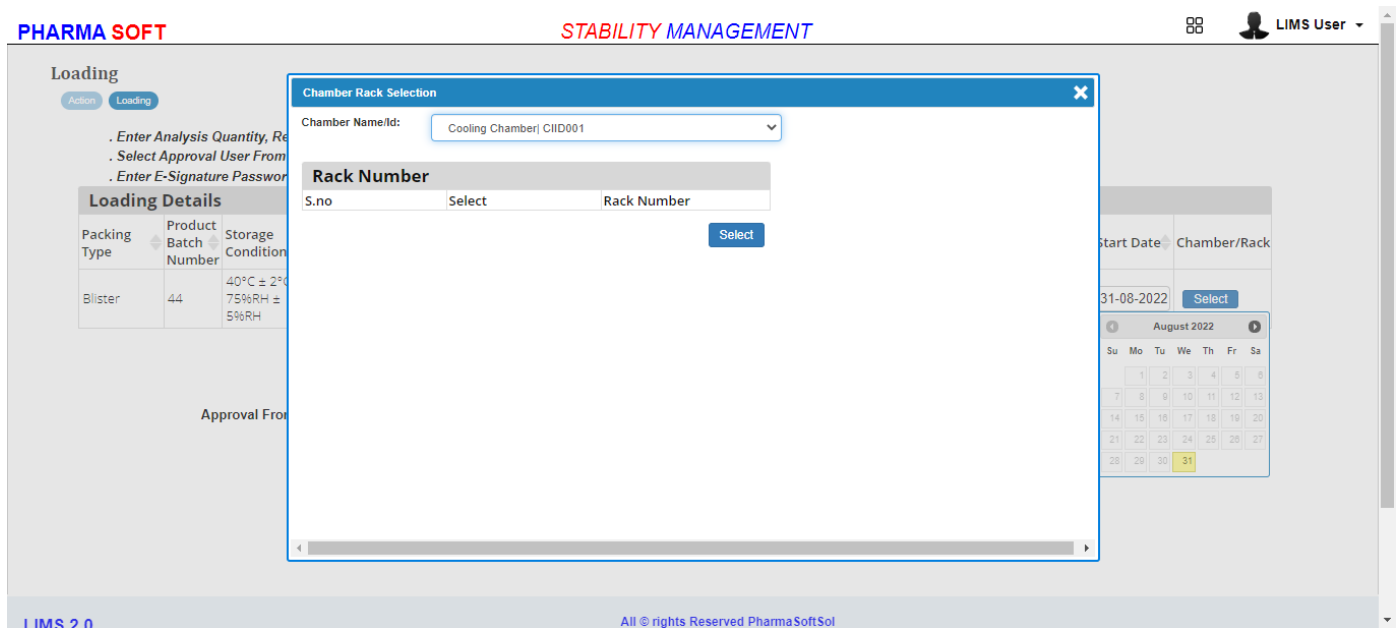

The screenshot shows the PHARMA SOFT STABILITY MANAGEMENT application. The 'Loading' section is active, displaying a table with 'Loading Details' including Packing Type (Blister), Product Batch Number (44), and Storage Condition (40°C ± 2°C, 75%RH ± 5%RH). A 'Chamber Rack Selection' dialog box is open, showing a 'Chamber Name/Id' dropdown set to 'Cooling Chamber| CIID001' and a 'Rack Number' table with columns 'S.no' and 'Rack Number'. A 'Select' button is visible in the dialog. In the background, a calendar for August 2022 is visible, with the date 31-08-2022 selected.

	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

13.Select Check Box

SSNo : PSS-LIMS-SS-013

PHARMA SOFT **STABILITY MANAGEMENT** LIMS User

Loading

Action Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
 . Select Approval User From Click On Select and Click On Finish.
 . Enter E-Signature Password and Click On Submit.

Loading Details

Packing Type	Product Batch Number	Storage Conditions
Blister	44	40°C ± 2°C / 75%RH ± 5%RH

Approval From:

Chamber Rack Selection

Chamber Name/Id:

S.no	Select	Rack Number
1	<input checked="" type="checkbox"/>	Rack-1
2	<input type="checkbox"/>	Rack-2

Start Date:

August 2022

Su	Mo	Tu	We	Th	Fr	Sa
1	2	3	4	5	6	
7	8	9	10	11	12	13
14	15	16	17	18	19	20
21	22	23	24	25	26	27
28	29	30	31			

IMS 2.0 All © rights Reserved PharmaSoftSol

14.Click on Select

SSNo : PSS-LIMS-SS-014

PHARMA SOFT **STABILITY MANAGEMENT** LIMS User

Loading

Action Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.
 . Select Approval User From Click On Select and Click On Finish.
 . Enter E-Signature Password and Click On Submit.

Loading Details


Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	<input type="text" value="2"/>	<input type="text" value="3"/>	Gm	<input type="text" value="3"/>	<input type="text" value="2"/>	Gm	<input type="button" value="Evaluate"/>	<input type="text" value="31-08-2022"/>	<input type="button" value="Select"/>

Approval From*:

August 2022


Su	Mo	Tu	We	Th	Fr	Sa
1	2	3	4	5	6	
7	8	9	10	11	12	13
14	15	16	17	18	19	20
21	22	23	24	25	26	27
28	29	30	31			

IMS 2.0 All © rights Reserved PharmaSoftSol

 PHARMA SOFT SOL <small>COMMITTED TO PAPERLESS COMPLIANCE</small>	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

17. Select a Record

SSNo : PSS-LIMS-SS-017

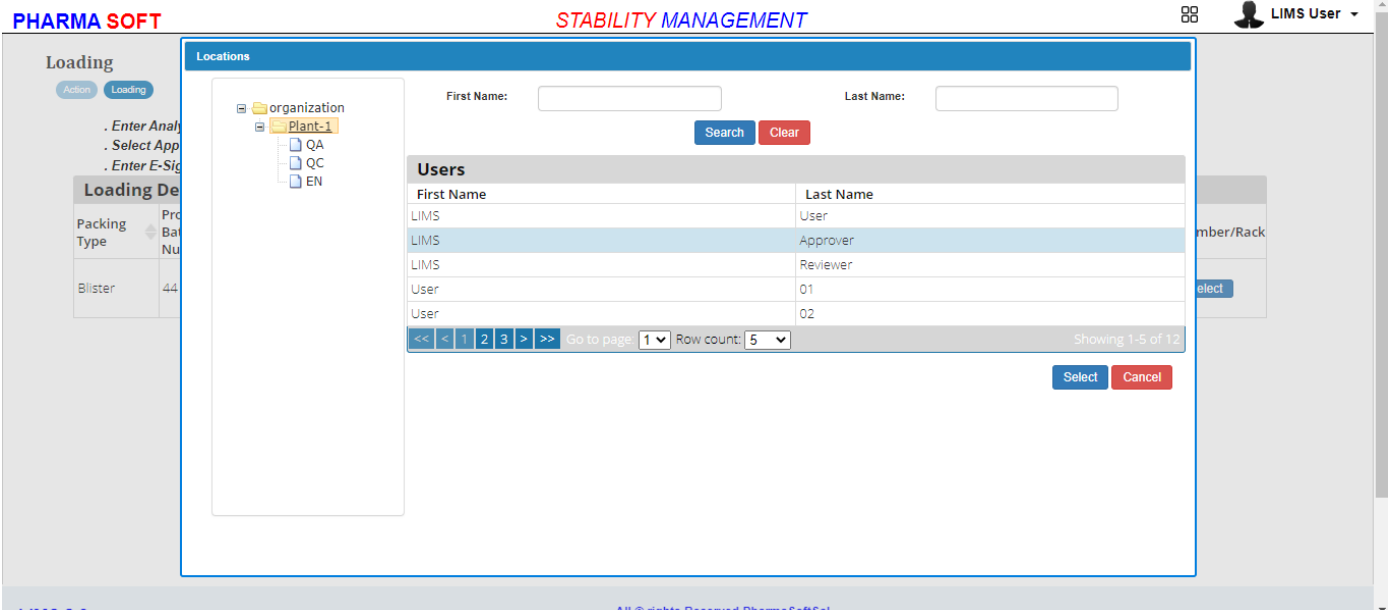

PHARMA SOFT **STABILITY MANAGEMENT** LIMS User

Loading

Action Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.

. Select Approval User From Click On Select and Click On Finish.

. Enter E-Signature Password and Click On Submit.

Loading Details


Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm	3	2	Gm	Evaluate	31-08-2022	Select

Approval From*: LIMS Approver Select

1 IMS 2 0 All © rights Reserved PharmaSoftSol

18. Click on Select

SSNo : PSS-LIMS-SS-018

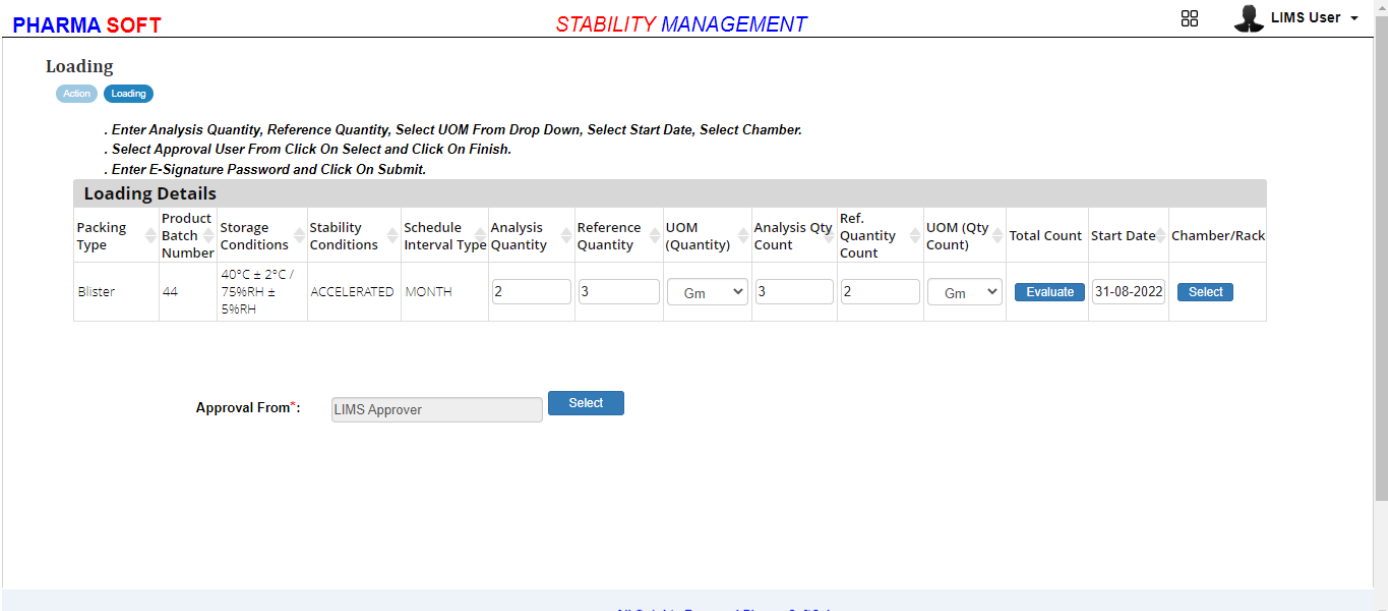

PHARMA SOFT **STABILITY MANAGEMENT** LIMS User

Loading

Action Loading

. Enter Analysis Quantity, Reference Quantity, Select UOM From Drop Down, Select Start Date, Select Chamber.

. Select Approval User From Click On Select and Click On Finish.


. Enter E-Signature Password and Click On Submit.

Loading Details

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm	3	2	Gm	Evaluate	31-08-2022	Select

Approval From*: LIMS Approver Select

1 IMS 2 0 All © rights Reserved PharmaSoftSol

	Test Case Name	Create Loading
	Test Case Number	LIMS-SM-011
	Actual Result	Pass

21.Click on Submit

SSNo : PSS-LIMS-SS-021

PHARMA SOFT **STABILITY MANAGEMENT** ☰ LIMS User

Enter E-Signature Password and Click on Submit

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm	3	2	Gm	Evaluate	31-08-2022	Select

Approval From*: LIMS Approver

Validate E-Signature

E-Signature Password*:

[Saving...](#)

[Submit](#) [Cancel](#)

[Previous](#) [Finish](#)

LIMS 2.0 All © rights Reserved PharmaSoftSol

22.Click on OK button

SSNo : PSS-LIMS-SS-022

PHARMA SOFT **STABILITY MANAGEMENT** ☰ LIMS User

Enter E-Signature Password and Click on Submit

Packing Type	Product Batch Number	Storage Conditions	Stability Conditions	Schedule Interval Type	Analysis Quantity	Reference Quantity	UOM (Quantity)	Analysis Qty Count	Ref. Quantity Count	UOM (Qty Count)	Total Count	Start Date	Chamber/Rack
Blister	44	40°C ± 2°C / 75%RH ± 5%RH	ACCELERATED	MONTH	2	3	Gm	3	2	Gm	Evaluate	31-08-2022	Select

Approval From*: LIMS Approver

Success

Submitted successfully.

[OK](#)

[Previous](#) [Finish](#)

LIMS 2.0 All © rights Reserved PharmaSoftSol

